

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	FORMIN		Page 1 of 71
BATCH SIZE	500,000 Tablets	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/014/00		

Generic Name:	Metformin Hydrochloride Tablets BP 500mg			BMR Revision No: 03	Approvals	Start date
Product Pack No:(√/ mark as appropriate)	FMN1TB	FMN2TJ	FMN3TJ	QA Head		
Pack Size: (√/ mark as appropriate)	10x10 Blisters	500's	1000's	Company/Prod. Pharmacist		
Reference SOP Number:	QAD/022/01			Production Manager		

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.***	Unit Quantity (mg/tablet)	Overage (mg/tablet)	Total* Quantity (mg/tablet)	Total Batch Quantity (Kg)	Lot	Quantity Per Lot	Tare Weight	Gross Weight	Net Weight	Scale ID.	Weighed by	Checked by**	Received by
1	Metformin Hydrochloride BP MET1HN		500.00	0.00	500.00	250.00	1	125.00kg							
							2	125.00kg							
2	Povidone K-30 BP PVP2XT		26.485	0.00	26.485	13.24	1	6.62kg							
							2	6.62kg							
3	Microcrystalline Cellulose BP MCC2XN		7.054	0.00	7.054	3.53	1	1.765Kg							
							2	1.765Kg							

PREPARED BY: Sign & Date:	REVIEWED BY: Sign & Date:	APPROVED BY: Sign & Date:	AUTHORIZED BY : Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	FORMIN		Page 2 of 71
BATCH SIZE	500,000 Tablets	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/014/00		

Generic Name:	Metformin Hydrochloride Tablets BP 500mg			BMR Revision No: 03	Approvals	Start date
Product Pack No:(√/ mark as appropriate)	FMN1TB	FMN2TJ	FMN3TJ	QA Head		
Pack Size: (√/ mark as appropriate)	10x10 Blisters	500's	1000's	Company/Prod. Pharmacist		
Reference SOP Number:	QAD/022/01			Production Manager		

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.***	Unit Quantity (mg/tablet)	Overage (mg/tablet)	Total* Quantity (mg/tablet)	Total Batch Quantity (Kg)	Lot	Quantity Per Lot	Tare Weight	Gross Weight	Net Weight	Scale ID.	Weighed by	Checked by**	Received by
4	Croscarmellose Sodium BP CMC2SN		10.582	0.00	10.582	5.29	1	5.29Kg							
5	Magnesium Stearate BP MAG2SN		3.527	0.00	3.527	1.764	1	1.764kg							
6	Purified Talc BP TAL2XN		2.351	0.00	2.351	1.176	1	1.176Kg							
TOTAL			550.00	0.00	550.00	0.00				Store-in charge	Dispensing Supervisor		QA Officer		
Weighing Started (HH:MM)			Date:		Weighing Stopped (HH:MM)			Date:		(Sign & Date)	(Sign & Date)		(Sign & Date)		
* Includes conversion to base drug, with overage and corrected for Potency and/or Moisture where required. ** Check weight, A.R. No .and verify presence of QC approved label on Raw Material container. *** AR No. Checked/Approved by QCM _____ (Sign & Date)															

PREPARED BY:	REVIEWED BY:	APPROVED BY:	AUTHORIZED BY :
Sign & Date:	Sign & Date:	Sign & Date:	Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	FORMIN		Page 3 of 71
BATCH SIZE	500,000 Tablets	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/014/00		

Generic Name:	Metformin Hydrochloride Tablets BP 500mg		BMR Revision No: 03	Approvals	Start date
Product Pack No:(√/ mark as appropriate)	FMN1TB	FMN2TJ	FMN3TJ	QA Head	
Pack Size: (√/ mark as appropriate)	10x10 Blisters	500's	1000's	Company/Prod. Pharmacist	
Reference SOP Number:	QAD/022/01			Production Manager	

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.***	Unit Quantity (mg/tablet)	Overage (mg/tablet)	Total* Quantity (mg/tablet)	Total Batch Quantity (Kg)	Lot	Quantity Per Lot	Tare Weight	Gross Weight	Net Weight	Scale ID.	Weighed by	Checked by**	Received by
7	Hydroxypropylmethyl Cellulose BP HPC2X1		4.00	0.00	4.00	2.00	1	500.00g							
							2	500.00g							
							3	500.00g							
							4	500.00g							
8	Propylene Glycol USP MPG2XN		0.53	0.00	0.53	0.265	1	66.25g							
							2	66.25g							
							3	66.25g							
							4	66.25g							

PREPARED BY: Sign & Date:	REVIEWED BY: Sign & Date:	APPROVED BY: Sign & Date:	AUTHORIZED BY : Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	FORMIN		Page 4 of 71
BATCH SIZE	500,000 Tablets	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/014/00		

Generic Name:	Metformin Hydrochloride Tablets BP 500mg		BMR Revision No: 03	Approvals	Start date
Product Pack No:(√/ mark as appropriate)	FMN1TB	FMN2TJ	FMN3TJ	QA Head	
Pack Size: (√/ mark as appropriate)	10x10 Blisters	500's	1000's	Company/Prod. Pharmacist	
Reference SOP Number:	QAD/022/01			Production Manager	

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.***	Unit Quantity (mg/tablet)	Overage (mg/tablet)	Total* Quantity (mg/tablet)	Total Batch Quantity (Kg)	Lot	Quantity Per Lot	Tare Weight	Gross Weight	Net Weight	Scale ID.	Weighed by	Checked by**	Received by
9	P.E.G. 6000 BP PEG2XN		0.60	0.00	0.60	0.30	1	75.00g							
							2	75.00g							
							3	75.00g							
							4	75.00g							
10	Isopropyl Alcohol BP IPA2XN		90.00	0.00	90.00	45.00	1	11.25Kg							
							2	11.25Kg							
							3	11.25Kg							
							4	11.25Kg							
11	Methylene Chloride MYC2XN		156.00	0.00	156.00	78.00	1	19.50Kg							
							2	19.50Kg							
							3	19.50Kg							
							4	19.50Kg							
12	Titanium Dioxide BP TIO3XN		1.14	0.00	1.14	0.57	1	142.50g							
							2	142.50g							
							3	142.50g							
							4	142.50g							

PREPARED BY: Sign & Date:	REVIEWED BY: Sign & Date:	APPROVED BY: Sign & Date:	AUTHORIZED BY : Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	FORMIN		Page 5 of 71
BATCH SIZE	500,000 Tablets	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/014/00		

Generic Name:	Metformin Hydrochloride Tablets BP 500mg		BMR Revision No: 03	Approvals	Start date
Product Pack No: (✓/ mark as appropriate)	FMN1TB	FMN2TJ	FMN3TJ	QA Head	
Pack Size: (✓/ mark as appropriate)	10x10 Blisters	500's	1000's	Company/Prod. Pharmacist	
Reference SOP Number:	QAD/022/01			Production Manager	

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.***	Unit Quantity (mg/tablet)	Overage (mg/tablet)	Total* Quantity (mg/tablet)	Total Batch Quantity (Kg)	Lot	Quantity Per Lot	Tare Weight	Gross Weight	Net Weight	Scale ID.	Weighed by	Checked by**	Received by
13	Talcum BP TAL2XN		0.50	0.00	0.50	0.25	1	62.50g							
		2					62.50g								
		3					62.50g								
		4					62.50g								
14	Talcum BP (for Polishing) TAL2XN		0.04	0.00	0.04	0.02	1	5.00g							
		2					5.00g								
		3					5.00g								
		4					5.00g								
Total			252.81	0.00	252.81	126.41				Store-in charge		Dispensing SPV		QA Officer	
Weighing Started (HH:MM) Date:				Weighing Stopped (HH:MM) Date:				(Sign & Date)		(Sign & Date)		(Sign & Date)			

* Includes conversion to base drug, with overage and corrected for Potency and/or Moisture where required.

** Check weight, A.R. No. and verify presence of QC approved label on Raw Material container.

*** AR No. Checked/Approved by QCM _____
(Sign & Date)

PREPARED BY: Sign & Date:	REVIEWED BY: Sign & Date:	APPROVED BY: Sign & Date:	AUTHORIZED BY : Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 6 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
PRODUCT DETAILS				

Label Claim	:	Each Formin Tablet Contains: Metformin Hydrochloride Tablets BP 500mg Colour: White
Shelf Life	:	2 years
Storage Conditions	:	Store in a cool & dry place. Protect from direct light, heat and moisture.
Mfg. License. Nos.	:	NDA/MAL/HDP/1701
Brand Name /Product Name	:	Formin
Required Batch Size	:	500,000 Tablets
Market	:	Domestic/Local
Average weight (Uncoated)	:	550.0mg
Average weight Coated	:	556.81mg
Dissolution	:	Not Less than 70 % for each Tablet
Assay:	:	Percentage Limit (95.0 % to 105.0 % of the Label claim)

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 7 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
INDEX				

Sr. Nos.	TITLES	Page Nos.
1.	Raw Material Requisition/Dispensing Sheets	1-6
2.	Product Details	7
3.	Index	8
4.	Abbreviations	9
5.	Equipment List& Instruments List	10
6.	General Instruction	11
7.	Safety Precaution	11
8.	General SOP In Manufacturing Respectively	11
9.	Dispensing Line Clearance	12
10.	Mixing/Granulation	13-18
11.	Blending/Lubrication	19-23
12.	Tablet Compression	24-41
13.	Tablet Inspection & Sorting	42-44
14.	Film Coating	45-50
15.	Tablet Inspection & Sorting	51-53
16.	Packaging	54-61
17.	In-Process Quality Control Report For Secondary Packing	64-66
18.	Finished Product Reconciliation Sheet	67
19.	Batch Release Sheet	68
20.	Revision History	69

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 8 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
ABBREVIATIONS				

S/N	SHORT	IN FULL	S/N	SHORT	IN FULL
01	BMR	Batch Manufacturing Record	19	L.O.D.	Loss On Drying
02	BPR	Batch Packing Record	20	NMT	Not More Than
03	API	Active Pharmaceutical Ingredients	21	NLT	Not Less Than
04	LAF	Laminar Air Flow	22	QAO	Quality Assurance Officer
05	RLAF	Reverse Laminar Air Flow	23	SPV	Supervisor
06	SOP	Standard Operating Procedure	24	IPQA	In-Process Quality Assurance
07	B/L	Break Line	25	IPA	Isopropyl Alcohol
08	AR No.	Analytical Report Number	26	NA	Not Applicable
09	Mfg.	Manufacturing	27	QA	Quality Assurance
10	Exp.	Expiry	28	M/C	Machine
11	RMG	Rapid Mixer Granulator	29	RPM	Revolutions Per Minute
12	FBD	Fluidized Bed Dryer	30	PMF No.	Product Master Formula Number
13	B. size	Batch Size	31	PVC	Polyvinyl Chloride
14	HDPE	High Density Polyethylene	32	BOPP	Biaxially Oriented Polypropylene
15	SMFPQ	Semi-Finished Product Quarantine	33	QC	Quality Control
16	IPQC	In-Process Quality Control	34	S/N	Serial Number
17	FIFO	First In-First Out	35	FEFO	First Expiry-First Out
18	LHS	Left Hand Side	36	RHS	Right Hand Side

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 9 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
LIST OF EQUIPMENTS & INSTRUMENTS					

List of Equipment used in the Manufacturing of Product.

S/N	Names	Mark	S/N	Names	Mark
DISPENSING					
1	Weighing Balance B-46		IPQC LAB		
2	Weighing Balance G-28				
3	Weighing Balance B-59				
4	RLAF Booth LAF-01		1	Weighing Balance B-53	
5	RLAF Booth LAF-05		2	Weighing Balance B-54	
GRANULATION 1			3	Weighing Balance B-55	
1	Rapid Mixer Granulator T-47		4	Weighing Balance B-23	
2	Fluidized Bed Dryers T-48		5	Weighing Balance B-52	
3	Fluidized Bed Dryers T-67		6	IR Moisture Analyzer B-19	
4	Vibro Sifter T-77		7	IR Moisture Analyzer B-43	
5	Paste Kettle T-49		8	Friability Apparatus QD-075	
6	Comminuting Mill T-50		9	Friability Apparatus T-33	
7	Comminuting Mill T-14		10	Hardness Tester QD-081	
8	Weighing Balance B-60		11	Hardness Tester QD-073	
GRANULATION 2			12	Disintegration Apparatus QD-084	
1	Rapid Mixer Granulator T-27		13	Disintegration Apparatus QD-085	
2	Fluidized Bed Dryers T-24		14	Vernier Caliper QD-086	
3	Fluidized Bed Dryers T-26		15	Vernier Caliper QD-087	
4	Vibro Sifter T-78		16	Leak Test Apparatus QD-088	
5	Paste Kettle T-60		17	Friability Apparatus QD-077	
COMPRESSION			VISUAL INSPECTION AND SORTING		
6	Multimill T-44		1	Metal Detector	
7	Comminuting Mill T-01		2	De-duster	
8	Weighing Balance B-13		BLISTER PACKING		
COMPRESSION			1	Blister Machine P-03	
1	Compression Machine T-06		2	Blister Machine P-10	
2	Compression Machine T-15		3	Blister Machine P-20	
3	Compression Machine T-30		4	Blister Machine P-21	
4	Compression Machine T-35		5	Blister Machine P-27	
5	Compression Machine T-52		STRIP PACKING		
6	Compression Machine T-80		1	Strip Packing Machine	
BLENDING			OTHERS (SPECIFY)		
1	Double Cone Blender T-09				
2	Double Cone Blender T-17				
3	Double Cone Blender T-04				
4	Vibrosifter T-66				
5	Vibrosifter T-79				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 10 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	

GENERAL INSTRUCTIONS

1. Use approved Equipment only, mentioned in respective stages of the BMR.
2. All Equipment for use must be duly cleaned as per respective SOPs and approved by QA.
3. Each container, machinery or Equipment must bear Status indicating label, duly signed by the person responsible.
4. All Equipment used, must have appropriate safety measures.
5. Care should be taken to ensure that all manufacturing Operations are performed as per the process under supervision of competent technical staff, workers should wear nose masks, and clean gloves, gowning throughout the manufacturing and packing operations.
6. Ensure that all the ingredients have Release labels for intact containers and proper A. R. Nos. for loose containers.
7. Ensure the materials are dispensed and checked as per respective SOPs and FIFO / FEFO system is followed.
8. Dispensing of materials should be performed in the presence of three qualified people, 1 from stores, 1 from Q.A. and 1 from production.
9. Observe the atmospheric conditions mentioned at various steps, in the BMR.
10. Dispensing labels shall be affixed with BMR.

SAFETY PRECAUTIONS

1. Safety Precautions to be taken to avoid explosions in FBD during drying.
2. Safety Precautions must be taken during grinding, milling, granulating process.
3. Protect Respiratory Organs by using cloth masks and gloves during manufacturing process.
4. Store all granules and tablets in well closed container lined with double polybags.
5. Follow personal hygienic requirements.
6. To avoid any accident, check the earthing and continuity of electricity while operating the machine.
7. Check the safety mechanisms at the start of the operations, at the end of operations, after power failures and any breakdowns of the machine.
8. **Material Transfers:** -Transfer all the dispensed materials from stores to manufacturing areas. Counter-check and record identities and weights of all the ingredients received. Do not use any material which does not comply with weight or analytical report (AR) numbers.

Storage of Materials: -Store the material, in process material, semi-finished and finished material in designated area properly labeled mentioning status of the material and batch details.

NOTE: Ensure all the SOPs used for references in the BMR are of current versions.

GENERAL SOPs IN MANUFACTURING RESPECTIVELY.		
S. Nos.	Titles	SOP Nos.
1.	Procedure for Environmental Monitoring Temperature & Differential Pressure.	QAD/040/00
2.	Procedure for Area Cleaning and sanitization of Premises.	GEN/008/06
3.	In-Process Quality Control and Assurance	QAD/021/00

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 11 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
DISPENSING					

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by
				QA (Sign & Date)
1. Area cleaning is done as per SOP.				
2. Ensure the absence of batch documents, labels, materials or remnants of previous product or batch.				
3. Ensure availability of the BMR for the in-coming batch and Verify that ALL materials and facilities for new batch are available, labeled and identified.				
4. Gowning and de-gowning procedure is followed.				
5. Dispensing Room, Dispensing booth, Balances and dispensing tools are clean and have status as ready for use.				
6. Clean scoops, containers and poly bags are available for dispensing as per SOP.				
7. Manometer reading of dispensing booth RLAF is within the required limit. Actual Manometer reading:				
8. Balance Calibration carried out before starting dispensing.				
9. Verify the release status & retest validity period of each material.				
10. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40-65)% _____				
11. Ensure that the RLAF is started 15 minutes before start of the dispensing activities.				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by
				QA (Sign & Date)
1. All material from current product have been removed.				
2. All facilities have been cleaned and labeled.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Ensure that utensils and accessories from previous operations have been removed.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 12 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
MIXING/GRANULATION					

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
				Status
1. Area cleaning as per SOP.				
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.				
3. Gowning procedure is followed.				
4. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.				
5. Check the proper cleanliness of Vibro-sifter, Paste Kettle, RMG, FBD, Multi mill etc.				
6. Ensure that the Vibrosifter sieves, Multimill screen, FBD bags, Scoops etc. cleaned.				
7. Check and ensure that Balance Calibration & Verification records are updated.				
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40-65)% _____				
9. Check the Integrity/Intactness of FBD gaskets, Product Bowl sieves and Dutch Meshes/Distribution Plates, and the sifting and milling sieves.				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
				Status
1. All material from current product have been removed.				
2. Check the cleanliness of Supply and Return air Grills, Light Fixtures, Waste Bins, and drainage points.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Machine(s) has been checked for cleanliness.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 13 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
MIXING/GRANULATION				

There are two lots of material for this batch, for each perform the following:

A. DRY MIXING

- Sift the following materials through mesh #20 fitted to a vibrator sifter into a double polythene-lined drum.
 - 125.00kg of Metformin Hydrochloride BP**
 - 1.765kg of Microcrystalline cellulose BP**
- Transfer the above materials into the Rapid mixer / granulator and mix for 10 minutes at SLOW speed and then at FAST speed for 5 minutes.

Lot	RMG NUMBER:								
	Speed	Impeller Speed [RPM]	Chopper Speed [RPM]	Started	Completed	Actual Time Taken (min)	Done by	Verified by SPV	Confirmed by QA
1	Slow								
	Fast								
2	Slow								
	Fast								

B. WET MIXING

Prepare Binders as follows:

- Sift 6.62Kg of povidone K-30 with sieve #40 and collect into double polythene-lined bag.
- Make slurry of sifted povidone K-30 with 3.25litres of purified water at room temperature in stainless steel.
- Boil 17.5 litres of purified Water in stainless steel vessel.
- Add the hot boiling water to the slurry in step 4 above while vigorously stirring. A smooth uniform paste forms.
- Add the binder in step 6 to the dry blend of step 2 and mix at SLOW speed for 10 minutes then at FAST speed for 10 minutes. Add additional amount of purified water with continuous mixing (if required).
- Unload the content into Fluid Bed Drier bowl.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 14 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
MIXING/GRANULATION					

Lot	RMG NUMBER:												
	Binder Preparation				Massing/Mixing								
	Specified Cold Water Volume	Actual Cold Water Volume	Specified Hot Water Volume	Actual Hot Water Volume	Speed	Impeller Speed [RPM]	Chopper Speed [RPM]	Started	Completed	Actual Time Taken (min)	Done by	SPV	QA
1					Slow								
					Fast								
2					Slow								
					Fast								

9. Dry the granules from step 8 above at 60 °C for 45 minutes in Fluid Bed Drier then remove and turn around the granules with a stainless steel scoop. Repeat the above procedure then mill the granules in Comminuting Mill fitted with sieve size 10. Redry at 60 °C for 45 minutes.

FIRST DRYING

LOT	FBD No.	DRYING CONDITIONS			DRYING TIME			Done by	SPV	QA
		STD Inlet Temp. (°C)	Observed Drying Temp. (°C)	Outlet Temp. (°C)	Started	Completed	Actual Time Taken (min)			
1										
2										

SECOND DRYING

LOT	FBD No.	DRYING CONDITIONS			DRYING TIME			Done by	SPV	QA
		STD Inlet Temp. (°C)	Observed Drying Temp. (°C)	Outlet Temp. (°C)	Started	Completed	Actual Time Taken (min)			
1										
2										

FINAL DRYING

Specifications: LOD Limits: _____ %.

LOT	FBD NO.	DRYING CONDITIONS			DRYING TIME				FINAL MILLING				
		Specified Inlet Temp.(°C)	Observed Drying Temp.(°C)	Outlet Temp.(°C)	Started	Completed	Actual Time Taken (min)	Miller No.	Sieve Size [mm]	LOD (%)	Done by	SPV	QA
1													
2													

10. Let the granules cool to room temperature then mill with Milling machine fitted with sieve No. 2.5. Check L.O.D on IR Balance at 105 °C (Limit 1-4 % w/w).

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 15 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
MIXING/GRANULATION					

QA REPORT

11. Loss on Drying: Check the moisture content in the granules:

Sample: Take 20g of dried granules and operate the moisture balance as per SOP.

Specifications: LOD Limits: _____%.

Acceptance criteria: % LOD should comply with specification:

Test Time	Lot no.	Drying Temp.(°C)	Drying duration (min)	FBD no.	LOD (%)	Moisture Balance ID	Remarks	Done by QA (Sign & Date)
	1							
	2							

Attach the slip from the moisture analyzer.

Release procedure: Yes / No

If the above tests fall within the specified limits, QA fills out the results in the above table and the material is sent to the next processing stage.

Rejection procedure: Yes / No

If the material fails any of the above tests, QA fills out and labels the RMG bowl with "REJECTED / REWORK" label and the material is sent for reworking.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 16 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
MIXING/GRANULATION				

12. Determine the percentage yield of the dried granules.

Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
Total			
Theoretical Weight			

Operator/ grp. Leader (sign & date) _____ Section Supervisor (sign & date) _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 17 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
MIXING/GRANULATION					

YIELD RECONCILIATION				
Steps	Descriptions	Qty. (Kg)	Done by SPV (Sign & Date)	Verified by QA (Sign & Date)
A.	Theoretical Batch Size with respect to dispensed material			
B.	Actual Quantities of Dried Granules			
C.	Sample Quantities (L.O.D/In-Process checks)			
D.	Validation Samples			
E.	Rejects (if any)			
F.	Total Actual Yield (B+C+D)			
G.	Unaccountable Losses [A-(E+F)]			
$\% \text{ yield} = F \times 100 / A = \dots\dots\dots\%$				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):
CAUSE OF VARIATION: YES/NO.

REMARKS/BREAKDOWNS.....
.....
.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 18 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
BLENDING/LUBRICATION					

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. Area cleaning as per SOP.					
2. Check out and ensure that all the previous product containers, materials, and labels are removed from the manufacturing area.					
3. Gowning procedure is followed.					
4. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
5. Check the proper cleanliness of the Air Supply and Return grills.					
6. Ensure that the Vibrosifter sieves, Double Cone Blender, Scoops etc. are cleaned.					
7. Check and ensure that Balance Calibration & Verification records are updated.					
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40-65)% _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. All material from current product have been removed.					
2. All facilities have been cleaned and labeled.					
3. The area has been cleaned and labeled.					
4. All paperwork for the current batch has been completed.					
5. Reconciliation of current batch has been done.					
6. Machine(s) has been checked for cleanliness.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 19 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
BLENDING/LUBRICATION					

1. Equipment used for Sifting and Blending

Equipment Name	Equipment ID No.	Mark (Equipment used)
Double Cone Blender	T-04	
Double Cone Blender	T-09	
Double Cone Blender	T-17	
Vibrosifter	T-66	
Vibrosifter	T-79	

2. Sift using Vibrosifter with sieves indicated in the table below then transfer into Double Cone Blender.

- _____ Kg Dried Granules.
- 1.764 Kg of Magnesium Stearate BP
- 5.29 Kg of Croscarmellose Sodium BP
- 1.176 Kg of Purified Talc BP

3. Add the Recoveries, if any (indicate quantities and batch Numbers)

Batch No.				
Quantity				

S/N	Ingredient	Sift through mesh size	Qty added (Kg)	Added by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
1	Dried Granules	#12				
2	Croscarmellose Sodium BP	#40				
3	Purified Talc BP	#40				
4	Recoveries	#40				
5	Magnesium Stearate BP	#40				
Total Weight of dry ingredients added						

4. Transfer into Double Cone Blender ingredients 1 to 4 in the table above and mix for 15 minutes.

5. Add the Magnesium Stearate to the blend in 4 above and mix for 5 minutes and unload into the double polybag lined containers. Record total number of containers, individual and total net weight of each container and label with all relevant product details.

Start time	Completion time	Time taken	specifications	Deviation	Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
Blending (Without Magnesium Stearate)							
Lubrication (With Magnesium Stearate)							

NB: If recoveries are added, Calculate the new batch size:

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 20 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
BLENDING/LUBRICATION					

Intermediate Sample for moisture content and homogeneity of blended / lubricated granules.

6. Take samples of blended granules as per SOP for sampling blended / lubricated granules and send to QA lab along with this record.

Quantity sampled	Top.....g	Middle.....g	Bottom.....g
Sampled by QA			

QA REPORT					
Test	Specification	Moisture Balance ID	Top Sample	Middle Sample	Bottom Sample
LOD by IR moisture balance (%)	1.0-4.0				

The samples **Comply** ☐ **Do not comply** ☐ With the above specifications.

QA Officer: (sign)..... Date.....

Release procedure: Yes / No

If the above tests fall within the specified limits, QA fills out and labels the blender with "PASSED" label and the material is sent to the next processing stage.

Rejection procedure: Yes / No

If the material fails any of the above tests, QA fills out and labels the blender with "REJECTED / REWORK" label and the material is sent for reworking.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 21 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
BLENDING/LUBRICATION				

7. Determine the percentage yield of the blended granules.

Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
Total			
Theoretical Weight			

Operator/ grp. Leader (sign & date) _____ Section Supervisor (sign & date) _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 22 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
BLENDING/LUBRICATION					

YIELD RECONCILIATION				
Steps	Descriptions	Qty. (Kg)	Done by SPV (Sign & Date)	Verified by QA (Sign & Date)
A.	Theoretical Batch Size with respect to Dried Granules + Blending Materials			
B.	Actual Quantities of Lubricated Granules			
C.	Sample Quantities (L.O.D/In-Process checks)			
D.	Validation Samples			
E.	Rejects (if any)			
F.	Total Actual Yield (B+C+D)			
G.	Unaccountable Losses [A-(E+F)]			
$\% \text{ yield} = F \times 100 / A = \dots\dots\dots\%$				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):
CAUSE OF VARIATION:

REMARKS/BREAKDOWNS.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 23 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. Area cleaning as per SOP.					
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.					
3. Gowning procedure is followed.					
4. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
5. Check the proper cleanliness of the Air Supply and Return grills.					
6. Ensure that the Vibrosifter sieves, Double Cone Blender, Scoops etc. are cleaned.					
7. Check and ensure that Balance Calibration & Verification records are updated.					
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40-65)% _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. All material from current product have been removed.					
2. All facilities have been cleaned and labeled.					
3. The area has been cleaned and labeled.					
4. All paperwork for the current batch has been completed.					
5. Reconciliation of current batch has been done.					
6. Machine(s) has been checked for cleanliness.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 24 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

1. EQUIPMENT USED FOR COMPRESSION

Equipment Name	Equipment ID	Mark	Equipment Name	Equipment ID	Mark
Compression Machine	T-06		Compression Machine	T-80	
Compression Machine	T-15		Metal Detector	T-84	
Compression Machine	T-30		Metal Detector	T-85	
Compression Machine	T-35		Deduster	T-86	
Compression Machine	T-52		Deduster	T-87	

2. MACHINE SET UP

TYPE OF TOOLING :-D-TOOLING		B-TOOLING	
Punch size: (for Trade)	12.5mm	Punch shape: (for Trade).	
Tablet face appearance/ description (Trade): White round biconvex uncoated tablets plain on both sides.		Upper:	Shallow Concave, Plain
		Lower:	Shallow Concave, Plain
Punch size: (for UG)	12.5mm	Punch shape: (for UG).	
Tablet face appearance/ description (UG): White round biconvex uncoated tablets plain on both sides.		Upper:	Shallow Concave, Plain
		Lower:	Shallow Concave, Plain
Cleaned by:	Date:	Machine speed:	
Machine setting done by:		Date:	
Checked by (section supervisor) :		Date:	
Previous Product:		Batch No.:	

3. Tablet Weight Range Limits

Individual tablet wt. (Nominal):			574.00mg
% of Nominal weight		Tablet wt. (mg)	Symbol
Upper range	>105%	>577.5	> +T2 Alert
	105%	577.5	= +T2 Alert
	103%	566.5	= +T1 Good
Norm	100%	550.0	Very Good
Lower range	97%	533.5	= -T1 Good
	95%	522.5	= -T2 Alert
	<95%	<522.5	< -T2 Alert

4. COMPRESSION

Started	Date:	Time:	Completed	Date:	Time:
---------	-------	-------	-----------	-------	-------

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 25 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET COMPRESSION				

5. DETERMINE PERCENTAGE YIELD

Date:	Shift:	Operator/ Grp. Leader	Weight of tablets compressed
Total weight of compressed tablets			Kg

YIELD RECONCILIATION

Steps	Descriptions	Qty. (Kg)	Done by SPV (Sign & Date)	Verified by QA (Sign & Date)
A.	Theoretical Batch Size with respect to Blended Granules			
B.	Actual Quantities of Tablets Compressed			
C.	Sample Quantities (In-Process checks)			
D.	Validation Samples			
E.	Rejects (if any)			
F.	Total Actual Yield (B+C+D)			
G.	Unaccountable Losses [A-(E+F)]			
% Yield =F x100/A =.....%				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):

CAUSE OF VARIATION:

REMARKS/BREAKDOWNS.....

.....

.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 26 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET COMPRESSION				

6. INITIAL COMPRESSION MACHINE SETTING INSTRUCTIONS:

NOTE: Collect the tablets of first three rounds and discard as non-recoverable residues.

1. Run feeder for one minute.
2. Set the weight of tablets as per limits given above by turning weight dossier anticlockwise for increasing weight and clockwise for decreasing weight & discard the tablets.
3. Check the weights of 20 tablets in group. Set weights by turning the weight dossier anticlockwise for increasing the weight and clockwise for decreasing the weight. Check total weights of twenty tablets three times & record. Calculate Average weight. Check the individual weights of all punches and record.
4. Initially check thickness of all stations tablets each side and during process.
5. Initially check hardness of all stations tablets each side and during process.
6. Check friability, disintegration time and individual weight variation of tablet.
7. Collect fresh samples after complete setting of parameters and visually inspect for any defects.
8. If it is found satisfactory then proceed for further process and verify all in process parameters and record as per next table.

STANDARD PARAMETERS		
PARAMETERS	FREQUENCIES	QUANTITIES
Appearance	Initial	All Station Tabs.
	Every 60 Minutes	20 Tablets/Side
Marking/ Embossing	Initial & Every 60 Minutes	All Station Tab.
Group Weight	Initial & Every 60 Minutes	20 Tablets/Side
Individual weight	Initial	All Station Tabs.
	Every 60 Minutes	20 Tablets/Side
Disintegration Time	Initial & Every 60 Minutes	6 Tablets/Side
Uniformity	Every 60 Minutes	-----
Hardness	Initial	All Station Tabs.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 27 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET COMPRESSION				

	Every 60 Minutes	All Station Tabs.
Thickness	Initial	All Station Tabs.
	Every 60 Minutes	All Station Tabs.
Friability	Every 60 Minutes	20 Tabs/side
M/c RPM	-----	-----

SPECIFICATIONS:

S/N	TESTS	SPECIFICATION
1.	Physical Appearance:	White round biconvex uncoated tablets plain on both sides.
2.	Diameter Limits (mm)	12.5 ±0.1
3.	Thickness Limits (mm)	5.00 ±0.3
4.	Hardness Limits (Kg/cm ²)	NLT 3
5.	Tablets Standard Weight (mg)	550.0 ± 2%
6.	Weight Uniformity Limits	± 5% of Average Weight. NMT 2 tablets should be out of 5% range and and none of the tablets should be out of 10% range.
7.	Friability Limits	NMT 1.0% If obviously cleaved, cracked, chipped or broken tablets are present in the sample after tumbling, the sample fails the test.
8.	Disintegration time (Min)	NMT 15

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 28 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

7. INITIAL DIES AND PUNCHES CHECKING RECORD

STA. No.	UPPER PUNCH	LOWER PUNCH	DIES SUITABLE	STA. No.	UPPER PUNCH	LOWER PUNCH	DIES SUITABLE	STA. No.	UPPER PUNCH	LOWER PUNCH	DIES SUITABLE
01				21				41			
02				22				42			
03				23				43			
04				24				44			
05				25				45			
06				26				46			
07				27				47			
08				28				48			
09				29				49			
10				30				50			
11				31				51			
12				32				52			
13				33				53			
14				34				54			
15				35				55			
16				36				56			
17				37				57			
18				38				58			
19				39				59			
20				40				60			

PUNCH SET BY (Operator Name)

CHECKED BY (Production Supervisor)

VERIFIED BY(QA).....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 29 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

8. INITIAL INDIVIDUAL WEIGHT RECORDS

Station. Nos.	Individual weights	Station. Nos.	Individual weights	Station. Nos.	Individual weights	Station. Nos.	Individual weights	Station. Nos.	Individual weights
LHS									
01		12		23		34		45	
02		13		24		35		46	
03		14		25		36		47	
04		15		26		37		48	
05		16		27		38		49	
06		17		28		39		50	
07		18		29		40		51	
08		19		30		41		52	
09		20		31		42		53	
10		21		32		43		54	
11		22		33		44		55	
RHS									
01		12		23		34		45	
02		13		24		35		46	
03		14		25		36		47	
04		15		26		37		48	
05		16		27		38		49	
06		17		28		39		50	
07		18		29		40		51	
08		19		30		41		52	
09		20		31		42		53	
10		21		32		43		54	
11		22		33		44		55	
Done by Operator (Sign & Date)			Checked by Supervisor (Sign & Date)			Verified by QA (Sign & Date)			

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 30 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

9. TABLET COMPRESSION IN-PROCESS QUALITY CONTROL REPORT									
M/C No.:				Location				Speed	
(I) Tablet Weight In-Process Control									
Time HH:MM									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Physical Appearance									
Individual Weights	1								
	2								
	3								
	4								
	5								
	6								
	7								
	8								
	9								
	10								
	11								
	12								
	13								
	14								
	15								
	16								
	17								
	18								
	19								
	20								
Total wt. (mg)									
Wt. after Friability									
Avg wt. (mg)									
Range									
Max. wt. (mg)									
Min. Wt. (mg)									
(II) Max/min wt. plot									
Action >+T2	>577.5								
Alert+T2	577.5								
Good +T1	566.5								
V. Good Norm	550.0								
Good -T1	533.5								
Alert -T2	522.5								
Action <-T2	<522.5								
Friability	NMT 1%								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 31 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

(III) Other In-Process Control Tests									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Disintegration	NMT 15min								
Hardness	NLT 3 (kg/cm ²)								
Diameter (mm)	12.50 ± 0.1								
Thickness (mm)	5.00 ± 0.3								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 32 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

TABLET COMPRESSION IN-PROCESS QUALITY CONTROL REPORT									
M/C No.:				Location				Speed	
(I) Tablet Weight In-Process Control									
Time HH:MM									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Physical Appearance									
Individual Weights	1								
	2								
	3								
	4								
	5								
	6								
	7								
	8								
	9								
	10								
	11								
	12								
	13								
	14								
	15								
	16								
	17								
	18								
	19								
	20								
Total wt. (mg)									
Wt. after Friability									
Avg wt. (mg)									
Range									
Max. wt. (mg)									
Min. Wt. (mg)									
(II) Max/min wt. plot									
Action >+T2	>577.5								
Alert+T2	577.5								
Good +T1	566.5								
V. Good Norm	550.0								
Good -T1	533.5								
Alert -T2	522.5								
Action <-T2	<522.5								
Friability	NMT 1%								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 33 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

(III) Other In-Process Control Tests									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Disintegration	NMT 15min								
Hardness	NLT 3 (kg/cm ²)								
Diameter (mm)	12.50 ± 0.1								
Thickness (mm)	5.00 ± 0.3								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 34 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

TABLET COMPRESSION IN-PROCESS QUALITY CONTROL REPORT									
M/C No.:				Location				Speed	
(I) Tablet Weight In-Process Control									
Time HH:MM									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Physical Appearance									
Individual Weights	1								
	2								
	3								
	4								
	5								
	6								
	7								
	8								
	9								
	10								
	11								
	12								
	13								
	14								
	15								
	16								
	17								
	18								
	19								
	20								
Total wt. (mg)									
Wt. after Friability									
Avg wt. (mg)									
Range									
Max. wt. (mg)									
Min. Wt. (mg)									
(II) Max/min wt. plot									
Action >+T2	>577.5								
Alert+T2	577.5								
Good +T1	566.5								
V. Good Norm	550.0								
Good -T1	533.5								
Alert -T2	522.5								
Action <-T2	<522.5								
Friability	NMT 1%								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 35 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

(III) Other In-Process Control Tests									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Disintegration	NMT 15min								
Hardness	NLT 3 (kg/cm ²)								
Diameter (mm)	12.50 ± 0.1								
Thickness (mm)	5.00 ± 0.3								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 36 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

TABLET COMPRESSION IN-PROCESS QUALITY CONTROL REPORT									
M/C No.:				Location				Speed	
(I) Tablet Weight In-Process Control									
Time HH:MM									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Physical Appearance									
Individual Weights	1								
	2								
	3								
	4								
	5								
	6								
	7								
	8								
	9								
	10								
	11								
	12								
	13								
	14								
	15								
	16								
	17								
	18								
	19								
	20								
Total wt. (mg)									
Wt. after Friability									
Avg wt. (mg)									
Range									
Max. wt. (mg)									
Min. Wt. (mg)									
(II) Max/min wt. plot									
Action >+T2	>577.5								
Alert+T2	577.5								
Good +T1	566.5								
V. Good Norm	550.0								
Good -T1	533.5								
Alert -T2	522.5								
Action <-T2	<522.5								
Friability	NMT 1%								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 37 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

(III) Other In-Process Control Tests									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Disintegration	NMT 15min								
Hardness	NLT 3 (kg/cm ²)								
Diameter (mm)	12.50 ± 0.1								
Thickness (mm)	5.00 ± 0.3								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 38 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

TABLET COMPRESSION IN-PROCESS QUALITY CONTROL REPORT									
M/C No.:				Location				Speed	
(I) Tablet Weight In-Process Control									
Time HH:MM									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Physical Appearance									
Individual Weights	1								
	2								
	3								
	4								
	5								
	6								
	7								
	8								
	9								
	10								
	11								
	12								
	13								
	14								
	15								
	16								
	17								
	18								
	19								
	20								
Total wt. (mg)									
Wt. after Friability									
Avg wt. (mg)									
Range									
Max. wt. (mg)									
Min. Wt. (mg)									
(II) Max/min wt. plot									
Action >+T2	>577.5								
Alert+T2	577.5								
Good +T1	566.5								
V. Good Norm	550.0								
Good -T1	533.5								
Alert -T2	522.5								
Action <-T2	<522.5								
Friability	NMT 1%								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 39 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET COMPRESSION					

(III) Other In-Process Control Tests									
Done by:		Operator:		QA:		Operator:		QA:	
Tablet discharge		Left	Right	Left	Right	Left	Right	Left	Right
Disintegration	NMT 15min								
Hardness	NLT 3 (kg/cm ²)								
Diameter (mm)	12.50 ± 0.1								
Thickness (mm)	5.00 ± 0.3								

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 40 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET COMPRESSION				

10. BULK TRANSFER / RECONCILIATION						
From: Tablet compression Section						
To: Semi-Finished product quarantine store (SMFPQ)						

Please receive the following tablets:

Date	Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)	Delivered by Operator	Received by SMFPQ Supervisor
	1					
	2					
	3					
	4					
	5					
	6					
	7					
	8					
	9					
	10					
Total						
Total number of tablets:				Average weight of tablets:		

Section supervisor (sign & date) _____

Remarks

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 41 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET INSPECTION & SORTING				

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date) Status
1. Area cleaning as per SOP.				
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.				
3. Gowning procedure is followed.				
4. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.				
5. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.				
6. Ensure that the Vibrosifter sieves, Metal Detector, Deduster, Scoops etc. are cleaned.				
7. Check and ensure that Balance Calibration & Verification records are updated.				
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date) Status
1. All material from current product have been removed.				
2. All facilities have been cleaned and labeled.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Machine(s) has been checked for cleanliness.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 42 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
TABLET INSPECTION & SORTING					

INSPECTION AND SORTING PROCEDURE

1. Ensure that the inspection area and equipment are clean.
2. Ensure that all personnel are properly dressed in protective garments, gloves and masks.
3. Carry out a line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
4. Bring the tablets inside the sorting room.
5. Visually inspect the tablets for defects such as chipping, spots etc.
6. Collect all rejected tablets in properly labeled designated containers.
7. Seal inspected tablets in polythene lined drums and properly label the drums.
8. Complete the "Ending activities" section of the line clearance form.
9. Transfer inspected and sorted tablets to Packing area/Semi-finished product quarantine (SMFPQ) store.

VISUAL INSPECTION LOG

ACTIVITY	FROM	TO

RECONCILIATION AFTER INSPECTION AND SORTING:

1	Weight of tablets received for inspection and sorting _____ Kg
2	Weight of tablets received after inspection and sorting _____ Kg
3	Average weight of tablets _____ mg/tablet
4	Recyclables after inspection and sorting _____ Kg (recyclables are recorded in the BMR reconciliation sheet)
5	Total rejects of tablets after inspection and sorting _____ Kg (rejects are recorded in the BMR reconciliation sheet).
6	Total loss after inspection and sorting _____ Kg _____ %

Done by **group leader (sign & date):**

Checked by **SMFPQ supervisor (sign & date):**

Verified by **QA (sign & date):**

7	Names of person(s) involved in carrying out Inspection and Sorting		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

NOTE: Transfer inspected and sorted tablets to Semi Finished Product Quarantine (SMFPQ) Store

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 43 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET INSPECTION & SORTING				

QA IN-PROCESS QUALITY CONTROL REPORT FOR INSPECTION & SORTING
--

Product compressed from Compression Machine Number _____

Record the net content of the drums to be sorted, weigh and record the quantity of defects sorted from the respective drums in the table below.

Drum number	Net weight (kg)	Quantity of defects (kg)	Nature of defects
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
Total			
Percentage defects = (Total quantity of defects / total net weight) x 100% = _____%			

Done by (QA officer, sign & date): _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 44 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
FILM COATING					

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.	Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
			QA (Sign & Date)	Status
1. Area cleaning as per SOP.				
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.				
3. Gowning procedure is followed.				
4. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.				
5. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.				
6. Ensure that the coating Pans, the spray gun assembly, silicon tubes and the gun stand, hot air blower pipes and exhaust pipes are cleaned.				
7. Check and Ensure that Balance Calibration & Verification records are updated.				
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements	Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
			QA (Sign & Date)	Status
1. All material from current product have been removed.				
2. All facilities have been cleaned and labeled.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Machine(s) has been checked for cleanliness.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 45 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
FILM COATING					

1. BULK TRANSFER / RECONCILIATION						
From: SMFPQ/ Inspection & Sorting Section						
To: Coating Section						
Please receive the following tablets:						
Date	Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)	Delivered by Operator	Received by SMFPQ Supervisor
	1					
	2					
	3					
	4					
	5					
	6					
	7					
	8					
	9					
	10					
Total						
Total number of tablets:				Average weight of tablets:		

Section supervisor (sign & date) _____

Remarks

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 46 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
FILM COATING				

COATING PRECAUTIONS

1. Once started the film coating process should not be interrupted.
2. The coating solutions should be kept for not more than 24 hours.
3. Once fixed the position of the gun should not be disturbed.
4. Compressed air should be free of water droplets.
5. Compressed air pressure should be maintained throughout film coating process.
6. Use a spray gun with a nozzle diameter of 1.0mm.

PREPARATION OF COATING DISPERSION

There are 4 lots for this batch, for every lot perform the following.

1. Place in a clean S.S. mixing vessel 16.50kg of Methylene Chloride.
2. In a separate S.S vessel, pour 3.0kg of Methylene Chloride, add while stirring 75.00gm of PEG 6000 to it and continue stirring until it completely dissolves then add in Step 1. above.
3. Add 66.25gm of Propylene Glycol in Step 1 while stirring continuously to get homogeneous solution.
4. Place in a clean S.S mixing vessel 4.50kg of Isopropanol Alcohol, add while stirring 500.00gm of HPMC to it and continue stirring until a smooth slurry is formed then add in step 1 above.
5. In another separate container pour about 6.75kg of Isopropanol Alcohol, add to it a mixture of 142.50gm of Titanium Dioxide and 62.50gm of Purified Talc, continue stirring until a smooth slurry is formed. Pass the solution through a fine sieve and pour the solution into the container in Step 1 while stirring.

COATING EQUIPMENT SETTINGS

S/N	PARAMETER	SET VALUE	ACTUAL READING
1	Spray gun pressure	2.0-2.25kg / cm ³	
2	Storage tank pressure	0.5-1.0kg/cm ³	
3	Coating pan speed	10-14 RPM	
4	Hot air blower Temperature	80 ± 5 °C	
5	Tablet bed Temperature	40 ± 5 °C	
Done by Operator (Sign & Date):			
Checked by SPV (Sign & Date):			
Verified by QA (Sign & Date):			

COATING PROCEDURE

1. Load the tablets into the pan and switch on the hot air blower and adjust the temperature of the incoming air at 80 °C. Jog the tablets for 15 min. by slow rotation of the pan intermittently until the tablet bed temperature is about 40 - 45 °C.
2. Switch on the exhaust fan and start rotating the pan locating the spray gun at the center of the tablet bed and to about 1/3 the height of the top of the falling tablet bed.
3. Spray the tablets with the coating solution without stopping until the solution is exhausted. Do not rotate the pan after the coating solution is exhausted.
4. Cool the tablets to about 40 °C then sprinkle with Talcum powder while rotating for about 5 minutes.
5. Off load the tablets into double polythene lined drums, weigh and close tightly.

NOTES

1. The coating solution should be kept for not more than 24 hours.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 47 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
FILM COATING					

1. TABLET COATING IN-PROCESS QUALITY CONTROL REPORT					
M/C No.:		Location		Speed	
(I) Tablet Weight In-Process Control					
Lot Number		1		2	
Time HH:MM					
Done by:		Operator	QA	Operator	QA
Physical Appearance					
Individual Weights	1				
	2				
	3				
	4				
	5				
	6				
	7				
	8				
	9				
	10				
	11				
	12				
	13				
	14				
	15				
	16				
	17				
	18				
	19				
	20				
Total wt. (mg)					
Avg wt. (mg)					
Range					
Max. wt. (mg)					
Min. Wt. (mg)					
Disintegration	NMT 30 min				
Hardness	NLT 4 (kg/cm ²)				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 48 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
FILM COATING					

TABLET COATING IN-PROCESS QUALITY CONTROL REPORT					
M/C No.:		Location		Speed	
(I) Tablet Weight In-Process Control					
Lot Number		3		4	
Time HH:MM					
Done by:		Operator	QA	Operator	QA
Physical Appearance					
Individual Weights	1				
	2				
	3				
	4				
	5				
	6				
	7				
	8				
	9				
	10				
	11				
	12				
	13				
	14				
	15				
	16				
	17				
	18				
	19				
	20				
Total wt. (mg)					
Avg wt. (mg)					
Range					
Max. wt. (mg)					
Min. Wt. (mg)					
Disintegration	NMT 30 min				
Hardness	NLT 4 (kg/cm ²)				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 49 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
FILM COATING					

DETERMINE PERCENTAGE YIELD OF COATED TABLETS

Date	Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)	Delivered by Operator	Received by SMFPQ Supervisor
	1					
	2					
	3					
	4					
	5					
	6					
	7					
	8					
Total						
Total number of tablets:				Average weight of tablets:		

YIELD RECONCILATION

Steps	Descriptions	Qty. (Kg)	Done by SPV (Sign & Date)	Verified by QA (Sign & Date)
A	Weight of tablets from SMFPQ			
B	Weight of coating materials			
C.	Theoretical Batch Size with respect to tablets received from SMFPQ + coating materials (A+B)			
D.	Actual Quantities of Tablets Coated			
E.	Sample Quantities (In-Process checks)			
F.	Validation Samples			
G.	Rejects (if any)			
H.	Total Actual Yield (D+E+F)			
I.	Unaccountable Losses [C-(G+H)]			
% Yield =H x100/C =.....%				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):

CAUSE OF VARIATION:

REMARKS/BREAKDOWNS.....

.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 50 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET INSPECTION & SORTING				

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
9. Area cleaning as per SOP.					
10. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.					
11. Gowning procedure is followed.					
12. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
13. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.					
14. Ensure that the Vibrosifter sieves, Metal Detector, Deduster, Scoops etc. are cleaned.					
15. Check and ensure that Balance Calibration & Verification records are updated.					
16. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
7. All material from current product have been removed.					
8. All facilities have been cleaned and labeled.					
9. The area has been cleaned and labeled.					
10. All paperwork for the current batch has been completed.					
11. Reconciliation of current batch has been done.					
12. Machine(s) has been checked for cleanliness.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 51 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET INSPECTION & SORTING				

INSPECTION AND SORTING PROCEDURE

1. Ensure that the inspection area and equipment are clean.
2. Ensure that all personnel are properly dressed in protective garments, gloves and masks.
3. Carry out a line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
4. Bring the tablets inside the sorting room.
5. Visually inspect the tablets for defects such as chipping, spots etc.
6. Collect all rejected tablets in properly labeled designated containers.
7. Seal inspected tablets in polythene lined drums and properly label the drums.
8. Complete the "Ending activities" section of the line clearance form.
9. Transfer inspected and sorted tablets to Packing area/Semi-finished product quarantine (SMFPQ) store.

VISUAL INSPECTION LOG		
ACTIVITY	FROM	TO
RECONCILIATION AFTER INSPECTION AND SORTING:		
1	Weight of tablets received for inspection and sorting _____ Kg	
2	Weight of tablets received after inspection and sorting _____ Kg	
3	Average weight of tablets _____ mg/tablet	
4	Recyclables after inspection and sorting _____ Kg (recyclables are recorded in the BMR reconciliation sheet)	
5	Total rejects of tablets after inspection and sorting _____ Kg (rejects are recorded in the BMR reconciliation sheet).	
6	Total loss after inspection and sorting _____ Kg _____ %	
Done by group leader (sign & date):		
Checked by SMFPQ supervisor (sign & date):		
Verified by QA (sign & date):		

7	Names of person(s) involved in carrying out Inspection and Sorting		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

NOTE: Transfer inspected and sorted tablets to Semi Finished Product Quarantine (SMFPQ) Store.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 52 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
TABLET INSPECTION & SORTING				

QA IN-PROCESS QUALITY CONTROL REPORT FOR INSPECTION & SORTING
--

Product coated from Coating Machine Number(s) _____

Record the net content of the drums to be sorted, weigh and record the quantity of defects sorted from the respective drums in the table below.

Drum number	Net weight (kg)	Quantity of defects (kg)	Nature of defects
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
Total			
Percentage defects = (Total quantity of defects / total net weight) x 100% = _____%			

Done by (QA officer, sign & date): _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 53 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	

PACKAGING MATERIALS REQUISITION SHEET

Requested by Production Supervisor			Name:		Signature:			Date:			
Approved by Production Manager			Name:		Signature:			Date:			
Pack Size: (/ mark as appropriate)					10x10 Blisters	100 Blisters X 10 Tablets	500's	1000's			
S/N	Item Code	Item Description	AR. NO	Total Batch Quantity	Units	Issued by (PM Stores In-charge)		Received by (Packaging Supervisor)		Verified by (QA)	
						Sign	Date	Sign	Date	Sign	Date
FOR BLISTER PACK											
1	2FMN1TB	Formin Tablets printed Aluminium Foil 203x0.25mm Trade)			Kg						
2	3FMN1TB	Formin Tablets printed Aluminium Foil 203 x0.25mm (UG)			Kg						
3	PF206M	PVC Clear film 206 x 0.25 mm			Kg						
4	1FMN1TB	Formin Tablets 10x10 Blister Cartons (Trade)			Kg						
5	4FMN1TB	Formin Tablets 10x10 Blister Cartons (UG)			PCS						
6	STB1XT	150x100 Tablets Blister Shipping Carton			PCS						
7	BPTPE2	2" BOPP Printed Packing Tape			PCS						
8	FMNPIL	Formin Tablets Patient Information Leaflet			PCS						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 54 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
PACKAGING				

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. Area cleaning as per SOP.					
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.					
3. Gowning procedure is followed.					
4. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
5. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.					
6. Ensure that the Vibrosifter sieves, Metal Detector, Deduster, Scoops etc. are cleaned.					
7. Remove all the cartons, labels, strips, any material left over on the line and on the lower space of the conveyor belts.					
8. Put an authorized overprinted specimen of strip/blister/cartons/Label/Tube in the BMR for ready reference.					
9. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. All material from current product have been removed.					
2. All facilities have been cleaned and labeled.					
3. The area has been cleaned and labeled.					
4. All paperwork for the current batch has been completed.					
5. Reconciliation of current batch has been done.					
6. Stereos of previous product are submitted to production Supervisor.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 55 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
PACKAGING				

1. PACKING PROCESS & INSTRUCTIONS

BULK PACKING

1. Ensure that the packing room, packing equipment and packing table(s) are clean.
2. Carry out a Packing line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
3. Bring in the appropriate packaging material.
4. Check the labels for correctness and legibility of coding.
5. Bring in the tablets to be packed.
6. Count/Weigh the tablets.
7. Put the tablets in poly bags and seal.
8. Place the sealed poly bag in the container and firmly place the cap.
9. Re-check the batch no., mfg. date and exp. date on the labels.
10. Carefully label the container.
11. Pack the right number of containers in the appropriate corrugated shipping cartons.
12. Seal the shipping cartons with packing tape.
13. Label the shipping cartons appropriately.
14. Complete the packing records.
15. Complete the "Ending activities" section of the packing line clearance form.
16. Transfer packed goods to Finished Goods Quarantine Store and wait for QC report.

BLISTER PACKING

1. Ensure that the packing room, packing equipment and packing table(s) are clean.
2. Carry out a Packing line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
3. Bring in the appropriate packaging material.
4. Set and operate the blister packing machine as per SOP for operating blister packing machine.
5. Check the first code for correctness and legibility of coding.
6. Bring in the tablets to be blistered and complete the operation and transfer all those blisters to secondary packing area.
7. Pack 10 blisters per unit carton/ 100 blisters per unit dispenser.
8. Pack 150 unit cartons in the corrugated shipping cartons.
9. Pack 24 dispenser in the corrugated shipping cartons.
10. Seal the shipping carton with packing tape.
11. Label the shipping carton appropriately.
12. Complete the packing records.
13. Complete the "Ending activities" section of the packing line clearance form.
14. Transfer packed goods to Finished Goods Quarantine Store and wait for QC report.

Reference Documents:

- 1.SOP for status labelling. 2. SOP for line clearance. 3. SOP for Bulk packing. 4. SOP for Blister packing

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 56 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
PACKAGING					

2. TABLET BLISTERING

1. BLISTER MACHINE SETTING

Blister Machine number		Blister machine Make	
PVC Clear Film Size		Aluminium Foil Size	
Set Sealing Temperature		Observed Sealing Temperature	
Set Forming Temperature		Observed Forming Temperature	
Machine Speed		Blister Pack Size	
Machine Setting Done By (Sign & Date):			
Checked by Supervisor (Sign & Date):			
Verified by QA (Sign & Date):			

2. BLISTER PACKING

Started	Date:	Time:	Completed	Date:	Time:
---------	-------	-------	-----------	-------	-------

3. DETERMINE PERCENTAGE YIELD

Date:	Shift:	Operator/ Grp. Leader	Number of Blisters
Total Number of Blisters			

YIELD RECONCILIATION

Steps	Descriptions	Qty. (Nos)	Done by SPV (Sign & Date)	Verified by QA (Sign & Date)
A.	Theoretical Batch Size (Blisters) with respect to tablets received from sorting			
B.	Actual Quantities of Blisters			
C.	Sample Quantities (In-Process checks)			
D.	Validation Samples			
E.	Rejects (if any)			
F.	Total Actual Yield (B+C+D)			
G.	Unaccountable Losses [A-(E+F)]			
Percentage Yield = $F \times 100 / A = \dots\dots\dots\%$				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):

CAUSE OF VARIATION: YES/NO

REMARKS/BREAKDOWNS.....

.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 57 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
PACKAGING					

3. PACKING CODING CONTROL

1	Batch No:	Coding approved by Packing Supervisor (sign & date)	
2	Mfg. Date:		
3	Exp.Date:		
4	Item(s) to be coded: Labels, Unit Cartons, Shipper labels, Blisters, Others (specify); mark as appropriate		
5	Coding Stamp set by (operator): sign & date:		
	Coding Machine Number:		
6	Attach first coded samples and patient information leaflet (PIL) here:		
7	Coding approved by: QA Officer (Sign & Date):		
8	Names of person(s) applying labels and packing - BULK & BLISTER (mark as appropriate).		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)
9	Coding Reconciliation:		
	Name of item:		
	No. of item issued:		
	No. of item coded:		
	No. of item uncoded:		
	No. of items damaged / rejected during coding:		
	No. of items used:		
	No. of excess items, if any:		
	No. of item to be destroyed (with QA consent):		
	Reconciliation approved by: QA (Sign & Date):		

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 58 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
PACKAGING					

4. BULK TRANSFER AFTER INSPECTION & SORTING

From: Semi finished product quarantine store (SMFPQ)

To: Packing Department (Bulk / blisters/strips)

Please receive the following tablets for packing:

Date	Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)	Delivered by SMFPQ supervisor	Received by Operator/ Packing supervisor
	1					
	2					
	3					
	4					
	5					
	6					
	7					
	8					
	9					
	10					
Total						
Total number of tablets:				Average weight of tablets:		

Quantity packed (bulk /blisters): -----

(A) REJECTS AFTER PACKING _____ KG/G

(B) RE-WORKABLES AFTER PACKING _____ KG/G

Section supervisor (sign & date) _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 62 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
SECONDARY PACKAGING					

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
9. Area cleaning as per SOP.					
10. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.					
11. Gowning procedure is followed.					
12. Ensure that status board is displayed mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
13. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.					
14. Ensure that the coating Pans, the spray gun assembly, silicon tubes and the gun stand, hot air blower pipes and exhaust pipes are cleaned.					
15. Check and Ensure that Balance Calibration & Verification records are updated.					
16. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
7. All material from current product have been removed.					
8. All facilities have been cleaned and labeled.					
9. The area has been cleaned and labeled.					
10. All paperwork for the current batch has been completed.					
11. Reconciliation of current batch has been done.					
12. Machine(s) has been checked for cleanliness.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 63 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
SECONDARY PACKAGING					

INSPECTION AND SORTING OF BLISTERS AND STRIPS			
INSPECTION AND SORTING PROCEDURE			
1. All personnel involved in Inspection, Sorting and Packing shall be qualified annually with eye vision checks (Ocular Test), colour-blindness checks, and ability to distinguish known product rejects from acceptable products. 2. The medical examination of all Inspectors/Sorters and Packers should include the report of proper vision with or without glasses and Color blindness. 3. The production Inspectors/Sorters and Packers who feel uneasy or complain of headache are unfit for handling Inspection, Sorting and Packing of Blisters/Strips. 4. Inspectors/Sorters who are employed for checking strips and blisters or overprinted materials shall be replaced with other trained checkers alternately after every four hours or a break of at least 30 minutes shall be taken to counter eye fatigue and complacency. 5. Before the start of any batch, line clearance should be done as per SOP for line clearance. 6. Bring the strips/blisters to the secondary packing room. 7. Visually inspect/Check the strips / blisters / coded materials for all the possible defects. Some examples of types of rejection on packing line are Coding rejection, Due to text matter/hologram/smudging of ink over printed foil, Rejected strips due to cutting, Empty pocket rejection, Pocket puncture rejection, Sealing rejection, Overheated strips/blisters, Chipped / capped / broken tablets, Improper forming of blisters, Improper appearance / absence of hologram in the final blister strip among others. 8. After the Inspection and Sorting operation, the Supervisor shall intimate to the QA Officer online to verify the integrity of the sorted blisters/strips and approve packing to start. 9. Remove all defective blisters/strips, place them in a polybag or drum and label as recyclables. These should be de-blistered/de-stripped once sorting and inspection of blisters is done.			
Inspection and Sorting Line Clearance Done by (sign & date)			
Checked by Production Supervisor (sign & date):			
Verified by QA (sign & date):			
VISUAL INSPECTION LOG			
ACTIVITY	FROM	TO	
10	Names of person(s) involved in carrying out Inspection and Sorting		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 64 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
SECONDARY PACKAGING					

SECONDARY PACKING PROCEDURE			
1. Check the specimen attached to the BMR of overprinted carton/Monocartons, leaflet, and shipper label if it has been checked by the line supervisor and approved by the QA Officer online. 2. Transfer previously approved overprinted packing material from the printing room to the secondary packing room. 3. Check if displayed details on the carton/Monocartons/dispensers of the batch number, manufacturing and expiry date are as per BMR/BPR. Check carton for overprinting intactness. 4. Check if the information on leaflets/inserts is clear, legible and specific to the product to be packed. 5. Open the carton and insert leaflet with number of strips/blisters as specified in the BMR. Close the carton and arrange it into the shipper. 6. Once the shipper is filled with the number of cartons as specified in the BMR, insert a fully filled packing slip with name of the product & batch number, the name of the packing personnel and this should be checked and signed by the group leader. 7. Group leaders should check the shipper to ensure that it has the correct quantity of cartons/Monocartons/dispensers. If this is okay, group leader shall sign the packing slip and intimate to the Supervisor and QA Officer online to confirm the correctness of the quantity in the shipper and other packing requirements. The Supervisor and the QA officer shall then sign on the packing slip and authorize the group leader to seal the shipper with BOPP tape and attach the shipper label. Exact quantity of loose shipper is written on shipper label and sealed. 8. Once the entire batch has been completed, transfer to finished product quarantine store. 9. Any leftover printed cartons/Monocartons/dispensers/labels should be destroyed in presence of QA Officer online and the number recorded in the BMR. 10. Transfer leftover packaging material (that is not printed on) to Packaging Material Stores with a return note.			
Secondary Packing Line Clearance Done by (sign & date)			
Checked by Production Supervisor (sign & date):			
Verified by QA (sign & date):			
RECONCILIATION AFTER PACKING			
Product Name:		Batch Number:	
Batch Size (Cartons/Jars):		Pack Size:	
No. of intact shippers packed		No. of cartons/jars in intact Shippers	
Loose Quantity Packed:		Total Quantity Packed (cartons/jars)	
Done by group leader (sign & date):			
Checked by SMFPQ supervisor (sign & date):			
Verified by QA (sign & date):			
11	Names of person(s) applying labels and packing - BULK & BLISTER (mark as appropriate).		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN		Page 65 of 71
BATCH SIZE		500,000 Tablets	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:	Shift :	
SECONDARY PACKAGING				

SHIPPER WEIGHT VERIFICATION

- After the packing operation is completed for the batch, **Weigh the packed shippers one by one** and record the weight of packed shipper in BPR and also on the Shipper Label and recheck the shipper if any discrepancy is found.

Average Weight of Empty Carton:					Average Weight of Packed Unit Carton:						
Average Weight of Empty Shipper:					Average Weight of Packed FULL Shipper:						
Acceptance Limit: _____ % of average weight of packed Shipper					Range: _____ to _____						
SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT
01		11		21		31		41		51	
02		12		22		32		42		52	
03		13		23		33		43		53	
04		14		24		34		44		54	
05		15		25		35		45		55	
06		16		26		36		46		56	
07		17		27		37		47		57	
08		18		28		38		48		58	
09		19		29		39		49		59	
10		20		30		40		50		60	

Remarks: All shippers are within Specifications YES/NO _____

Weighing Done by (Sign & Date):

Checked by Packing Supervisor (Sign & Date):

Approved by QA (Sign & Date):

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 69 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Shift :	
FINISHED PRODUCT RECONCILIATION SHEET					

Quantity Packed	Retention Sample Quantity	Actual Quantity Released	Quantity Delivered to FGS	Delivery Note Number	Delivery Date

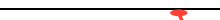
FINAL BATCH RECONCILIATION			
ITEM	Quantity (Gm/Kg).	Quantity* (No. of Units)	Percent =[QTYx100/(A)]
(A) Theoretical Yield =(Batch Size):			
(B) Yield - 1st Crop (Quantity Packed):			
(C) In-Process/QC Samples:			
(D) Re-workables =[(I)+(II)+ (III)]:			
(I) Powder/Granules*			
(II)Uncoated Tablets			
(III) Coated Tablets			
(E) Rejects =[(I)+(II)+(III)]:			
(I) Powder/Granules*			
(II) Uncoated Tablets			
(III) Coated Tablets			
(F) Total Yield =(B+C+D):			
(H) Shortage =[A-(E+F)]:			
Done by Supervisor (Sign & Date):		Verified by QAO (Stamp, Sign & Date):	

Reconciliation Approved By **Production Manager (Sign & Date):** _____

*Quantity (No. of units) is derived by converting this quantity using average compression weight.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		FORMIN		Page 70 of 71	
BATCH SIZE		500,000 Tablets		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/014/00	Date:		Actual Quantity :	
BATCH RELEASE SHEET					

The following Batch production records have been properly filled out, checked, dated and signed by authorized person(s). The tests are also within limits. The batch is hereby approved for release.		
Reviewed and approved by QA Head/Designee (Sign & Date):		
Records Reviewed (Mark as appropriate)	QA Head/Designee	Company /Prod. Pharmacist
1. Raw Material Requisition/Dispensing Sheet		
2. Mixing / Granulation Records		
3. Blending / Lubrication Records		
4. Compression/ In-process Control Records		
5. Process Line Clearance Records		
6. Records for Semi-finished Products (blended material, compressed tablets etc.)		
7. Inspection and Sorting Records		
8. Packing Line Clearance Records		
9. Packing Material Authorization Note		
10. Packing Coding Control Records		
11. Product Reconciliation Sheet		
12. IPQA Reports		
13. QC Analytical reports		
13.1. Finished Product Analytical Report		
13.2. Finished Product Analytical Data Sheet (Observations, Calculations, Printouts, References).		
I hereby certify that all the manufacturing stages of this batch of finished product have been carried out in full compliance with the Uganda National Drug Authority (NDA) GMP requirements and with the requirements of the product registration. I, therefore, AUTHORISE release of the batch.		
Company Pharmacist/Pharmacist In-Charge:		
_____	_____	_____
Name	Signature	Date

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		FORMIN	Page 71 of 71	
BATCH SIZE		500,000 Tablets		
MFR No.	KPI/MFR/014/00			

REVISION NO.	CHANGES INCORPORATED	REASON OF CHANGE	EFFECTIVE DATE
01	- KPI new logo to be incorporated - The provision for BMR preparation, Review, Approval and Authorization have been introduced on the Master Formula Page(s). - An Index has been introduced. - The Abbreviations used in the BMR have been defined in full. - List Of Equipment & Instruments used in product manufactured have been summarized. - Line Clearance for Dispensing previously missing has been introduced. - Documentation of First Drying, Second Drying, and Final Drying have been split and captured independently, and Milling activities including Sieves used plus LOD data are also being documented at granulation. Identity of the FBD used to dry a particular Lot has been captured. - In-process Quality Control and Assurance (IPQC/IPQA) reports for both Production and QA at all stages (as appropriate) have been incorporated. - Initial Compression Machine Setting instructions, Initial Dies and Punches Checking Record, Initial Individual Weight Records as per number of Punches/Stations of the Compression machine have been introduced/Captured.	To improve on clarity of the BMR. To match up with the current practices.	
02	- To eliminate duplication in the BMR header, introduce page breaks for different process steps and include MFR numbers in the BMR. - To include the Line Clearance for Inspection and Sorting of Blisters/Strips and the Secondary Packing Procedure.	To improve coherence of the BMR To ensure maximum elimination of blister/Strip defects before packing.	
03	- Batch Release sheet revised to ensure that batch release is Authorized by the Company Pharmacist/Pharmacist In-Charge. - Shipper weight verification record introduced. - Packaging Material Requisition Sheet incorporated in the BMR.	To comply with regulatory requirements for batch release. To eliminate the risk of having less or more product quantity in a shipper	