

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	KAM AMOXY		Page 1 of 45
BATCH SIZE	500,000 Capsules	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/033/00		

Generic Name:	Amoxicillin Capsules BP 250mg	BMR Revision No: 03	Approvals	Start date
Product Pack No:(√/ mark as appropriate)	KMX1CB	KMX2CB	QA Head	
Pack Size: (√/ mark as appropriate)	10x10 Blisters	1000's	Company/Prod. Pharmacist	
Reference SOP Number:	QAD/022/01		Production Manager	

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.****	Unit Quantity (mg/Capsule)	Overage (mg/Capsule)	Total Quantity** (mg/Capsule)	Total Batch Quantity (kg)	Lot	Quantity Per Lot	Tare Weight (Kg/g)	Gross Weight (Kg/g)	Net Weight (Kg/g)	Scale No.*	Weighed by	Checked by***	Received by
1	Amoxicillin Trihydrate BP Compacted AMX1TC		288.00	2.00	290.00	145.0kg	1	145.00kg							
2	Colloidal Silicone Dioxide BP (Aerosil 200) AER2XN		2.00	0	2.00	1.0kg	1	1.00kg							
3	Lactose Monohydrate BP LAC2MN		21.00	0	21.00	10.5kg	1	10.5kg							
4	Magnesium Stearate BP MAG2SN		7.00	0	7.00	3.5kg	1	3.5kg							
TOTAL			318.00	2.00	320.00	160.000kg									

PREPARED BY: Sign & Date:	REVIEWED BY: Sign & Date:	APPROVED BY: Sign & Date:	AUTHORIZED BY : Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

QAD/022/F01/01

	Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT	KAM AMOXY		Page 2 of 45
BATCH SIZE	500,000 Capsules	Mfg. Date :	
BATCH No.		Exp. Date :	
Effective Date:	MFR No.: KPI/MFR/033/00		

Generic Name:	Amoxicillin Capsules BP 250mg	BMR Revision No: 03	Approvals	Start date
Product Pack No:(√/ mark as appropriate)	KMX1CB	KMX2CB	QA Head	
Pack Size: (√/ mark as appropriate)	10x10 Blisters	1000's	Company/Prod. Pharmacist	
Reference SOP Number:	QAD/022/01		Production Manager	

RAW MATERIAL REQUISITION/DISPENSING SHEETS

Sr.No.	Description	A.R. No.****	Unit Quantity (mg/Capsule)	Overage (mg/Capsule)	Total Quantity** (mg/Capsule)	Total Batch Quantity (kg)	Lot	Quantity Per Lot	Tare Weight (Kg/g)	Gross Weight (Kg/g)	Net Weight (Kg/g)	Scale No.*	Weighed by	Checked by***	Received by
5	Empty Hard Gelatin Shell # 1 Maroon/Peach Printed KAMAMOXY 250 (body), KPI (cap) GMP5HN					510,000	1	510,000							

Weighing Started (HH:MM)	Date:	Weighing Stopped (HH:MM)	Date:	Store In-charge	Dispensing SPV	QA Officer
				(Sign & Date)	(Sign & Date)	(Sign & Date)

** Includes conversion to base drug, with overage and corrected for Potency and/or Moisture where required.

*** Check weight, A.R. No .and verify presence of QC approved label on Raw Material container.

**** AR No. Checked/Approved by QCM _____
(Sign & Date)

PREPARED BY: Sign & Date:	REVIEWED BY: Sign & Date:	APPROVED BY: Sign & Date:	AUTHORIZED BY : Sign & Date:
Name:	Name:	Name:	Name:
Designation:	Designation:	Designation:	Designation:

QAD/022/F01/01

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 3 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
PRODUCT DETAILS				

Label Claim	:	Each Kam Amoxy Capsule Contains: Amoxicillin Capsules BP 250mg Colour: Cap: Maroon printed on ‘KPI logo’/ UG & KPI logo in black. Body: Peach printed on ‘KAM AMOXY 250’ in black
Shelf Life	:	3 years
Storage Conditions	:	Store in a cool & dry place. Protect from direct light, heat and moisture.
Mfg. License Nos.	:	NDA/MAL/HDP/8089
Brand Name /Product Name	:	Kam Amoxy
Required Batch Size	:	500,000 Capsules
Market	:	Domestic/Local
Average weight (Net Content)	:	320.0mg
Dissolution	:	NA
Assay:	:	Percentage Limit (92.5 % to 110.0 % of the Label claim)

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 4 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
INDEX				

Sr. Nos.	TITLES	Page Nos.
1.	Raw Material Requisition/Dispensing Sheets	1-2
2.	Product Details	3
3.	Index	4
4.	Abbreviations	5
5.	Equipment List& Instruments List	5
6.	General Instruction	6
7.	Safety Precaution	6
8.	General SOP In Manufacturing Respectively	6
9.	Dispensing	7
10.	Blending/Lubrication	8-12
11.	Capsule Filling	13-23
12.	Capsule Inspection & Sorting	24-26
13.	Packaging	27-32
14.	In-Process Quality Control Report For Primary Packing	33-35
15.	In-Process Quality Control Report For Secondary Packing	36-38
16.	Finished Product Reconciliation Sheet	41
17.	Batch Release Sheet	42
18.	Revision History	43

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 5 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
ABBREVIATIONS					

S/N	SHORT	IN FULL	S/N	SHORT	IN FULL
01	BMR	Batch Manufacturing Record	19	L.O.D.	Loss On Drying
02	BPR	Batch Packing Record	20	NMT	Not More Than
03	API	Active Pharmaceutical Ingredients	21	NLT	Not Less Than
04	LAF	Laminar Air Flow	22	QAO	Quality Assurance Officer
05	RLAF	Reverse Laminar Air Flow	23	SPV	Supervisor
06	SOP	Standard Operating Procedure	24	IPQA	In-Process Quality Assurance
07	B/L	Break Line	25	IPA	Isopropyl Alcohol
08	AR No.	Analytical Report Number	26	NA	Not Applicable
09	Mfg.	Manufacturing	27	QA	Quality Assurance
10	Exp.	Expiry	28	M/C	Machine
11	RMG	Rapid Mixer Granulator	29	RPM	Revolutions Per Minute
12	FBD	Fluidized Bed Dryer	30	PMF No.	Product Master Formula Number
13	B. size	Batch Size	31	PVC	Polyvinyl Chloride
14	HDPE	High Density Polyethylene	32	BOPP	Biaxially Oriented Polypropylene
15	SMFPQ	Semi-Finished Product Quarantine	33	QC	Quality Control
16	IPQC	In-Process Quality Control	34	S/N	Serial Number
17	FIFO	First In-First Out	35	FEFO	First Expiry-First Out
18	LHS	Left Hand Side	36	RHS	Right Hand Side

LIST OF EQUIPMENTS & INSTRUMENTS

List of Equipment used in the Manufacturing of Product.

S/ N	Equipment Name	Equipment No.	Mark	S/ N	Equipment Name	Equipment No.	Mark
1	Multimill	PN-09		12	Leak Test Machine	QD-070	
2	Double Cone Blender	PN-06		13	LAF Booth		
3	Mechanical Sifter	PN-07		14	Weighing Scale		
4	Compression Machine	PN-35		15	Capsule Filling Machine	PN-37	
4	Tablet Dedusting Machine	PN-30		16	Capsule Mini Sorter	PN-44	
5	Tablet Dedusting Machine	PN-31		17	Capsule Polishing Machine	PN-40	
6	Blister Packing Machine	P-04		18	Capsule Loading Machine	PN-42	
7	Hardness tester	QD-082		19	Strip Packing Machine		
8	Friabilator	QD-080		20	Vernier Caliper		
9	Disintegration Tester Apparatus	QD-091		OTHERS (SPECIFY)			
10	Electronic Moisture Balance	B-20					
11	Electronic Balance	B-33					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 6 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	

GENERAL INSTRUCTIONS

1. Use approved Equipment only, mentioned in respective stages of the BMR.
2. All Equipment for use must be duly cleaned as per respective SOPs and approved by QA.
3. Each container, machinery or Equipment must bear Status indicating label, duly signed by the person responsible.
4. All Equipment used, must have appropriate safety measures.
5. Care should be taken to ensure that all manufacturing Operations are performed as per the process under supervision of competent technical staff, workers should wear nose masks, and clean gloves, gowning throughout the manufacturing and packing operations.
6. Ensure that all the ingredients have Release labels for intact containers and proper A. R. Nos. for loose containers.
7. Ensure the materials are dispensed and checked as per respective SOPs and FIFO / FEFO system is followed.
8. Dispensing materials should be performed in the presence of three qualified people, 1 from stores, 1 from Q.A. and 1 from production.
9. Observe the atmospheric conditions mentioned at various steps, in the BMR.
10. Dispensing labels shall be affixed with BMR.

SAFETY PRECAUTIONS

1. Safety Precautions to be taken to avoid explosions in FBD during drying.
2. Safety Precautions must be taken during grinding, milling, granulating process.
3. Protect Respiratory Organs by using cloth masks and gloves during manufacturing process.
4. Store all granules and tablets in well closed container lined with double polybags.
5. Follow personal hygienic requirements.
6. To avoid any accident, check the earthing and continuity of electricity while operating the machine.
7. Check the safety mechanisms at the start of the operations, at the end of operations, after power failures and any breakdowns of the machine.
8. **Material Transfers:** -Transfer all the dispensed materials from stores to manufacturing areas. Counter-check and record identities and weights of all the ingredients received. Do not use any material which does not comply with weight or analytical report (AR) numbers.

Storage of Materials: -Store the material, in process material, semi-finished and finished material in designated area properly labeled mentioning status of the material and batch details.

NOTE: Ensure all the SOPs used for references in the BMR are of current versions.

GENERAL SOPs IN MANUFACTURING RESPECTIVELY.		
S. Nos.	Titles	SOP Nos.
1.	Procedure for Environmental Monitoring Temperature & Differential Pressure.	QAD/040/00
2.	Procedure for Area Cleaning and sanitization of Premises.	GEN/008/06
3.	In-Process Quality Control and Assurance	QAD/021/00

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 7 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
DISPENSING					

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. Area cleaning is done as per SOP.					
2. Ensure the absence of batch documents, labels, materials or remnants of previous product or batch.					
3. Ensure availability of the BMR for the in-coming batch and Verify that ALL materials and facilities for new batch are available, labeled and identified.					
4. Gowning and de-gowning procedure is followed.					
5. Dispensing Room, Dispensing booth, Balances and dispensing tools are clean and have status as ready for use.					
6. Clean scoops, containers and poly bags are available for dispensing as per SOP.					
7. Manometer reading of dispensing booth RLAF is within the required limit. Actual Manometer reading:					
8. Balance Calibration carried out before starting dispensing.					
9. Verify the release status & retest validity period of each material.					
10. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40-65)% _____					
11. Ensure that the RLAF is started 15 minutes before start of the dispensing activities.					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. All material from current product have been removed.					
2. All facilities have been cleaned and labeled.					
3. The area has been cleaned and labeled.					
4. All paperwork for the current batch has been completed.					
5. Reconciliation of current batch has been done.					
6. Ensure that utensils and accessories from previous operations have been removed.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 8 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
BLENDING/LUBRICATION					

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
				Status
1. Area cleaning as per SOP.				
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.				
3. Gowning procedure is followed.				
4. Ensure that status board is displayed with mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.				
5. Check the proper cleanliness of the Air Supply and Return grills.				
6. Ensure that the Vibrosifter sieves, Double Cone Blender, Scoops etc. are cleaned.				
7. Check and ensure that Balance Calibration & Verification records are updated.				
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 - 65)% _____				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
				Status
1. All material from current product have been removed.				
2. All facilities have been cleaned and labeled.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Machine(s) has been checked for cleanliness.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 9 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
BLENDING/LUBRICATION					

1. Equipment used for Sifting and Blending

Equipment Name	Equipment ID No.	Mark (Equipment used)
Double Cone Blender	PN-06	
Multimill	PN-09	
Mechanical Sifter	PN-07	

2. Sift using mechanical Sifter with sieves indicated in the table below then transfer into Double Cone Blender.

- 145.00kg Amoxicillin Trihydrate Compacted
- 1.00kg Colloidal Silicon Dioxide BP (Aerosil 200)
- 10.5kg Lactose Monohydrate BP
- 3.5kg Magnesium Stearate BP
- _____ Powder from returned Capsules

3. Add the Recoveries, if any (indicate quantities and batch Numbers)

Batch No.				
Quantity				

4. Transfer into Double Cone Blender approximately one half the quantity of item (a), then add items (b), (c) and (d) and finally the other half of (a) and fill in the table below.

S/N	Ingredient	Sift through mesh size	Qty added (Kg)	Added by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
1	Amoxicillin Trihydrate Compacted	#20				
2	Colloidal Silicone Dioxide	#20				
3	Lactose Monohydrate	#20				
4	Magnesium Stearate BP	#20				
5	Powder from returned Capsules	#20				
Total Weight of dry ingredients added						

5. After all ingredients have been added to the double cone blender, Blend for 15 minutes and fill in the table below.

Start time	Completion time	Time taken	specifications	Deviation	Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)

NB: If recoveries are added, Calculate the new batch size:

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 10 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
BLENDING/LUBRICATION					

Intermediate Sample for moisture content and homogeneity of blended / lubricated granules.

6. Take samples of blended granules as per SOP for sampling blended / lubricated granules and send to QA lab along with this record.

Quantity sampled	Top.....g	Middle.....g	Bottom.....g
Sampled by QA			

QA REPORT					
Test	Specification	Moisture Balance ID	Top Sample	Middle Sample	Bottom Sample
LOD by IR moisture balance (%)	NMT 1.5				

The samples **Comply** ☐ **Do not comply** ☐ With the above specifications.

QA Officer: (sign)..... Date.....

Release procedure: Yes / No

If the above tests fall within the specified limits, QA fills out and labels the blender with "PASSED" label and the material is sent to the next processing stage.

Rejection procedure: Yes / No

If the material fails any of the above tests, QA fills out and labels the blender with "REJECTED / REWORK" label and the material is sent for reworking.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 11 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
BLENDING/LUBRICATION				

7. Determine the percentage yield of the blended granules.

Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
Total			
Theoretical Weight			

Operator/ grp. Leader (sign & date) _____ Section Supervisor (sign & date) _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 12 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
BLENDING/LUBRICATION					

YIELD RECONCILIATION				
Steps	Descriptions	Qty. (Kg)	Done by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
A.	Theoretical Batch Size with respect to Dispensed Materials			
B.	Actual Quantities of Lubricated Granules			
C.	Sample Quantities (L.O.D/In-Process checks)			
D.	Validation Samples			
E.	Rejects (if any)			
F.	Total Actual Yield (B+C+D)			
G.	Unaccountable Losses [A-(E+F)]			
Percentage yield = $F \times 100 / A = \dots\dots\dots\%$				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):
CAUSE OF VARIATION:

REMARKS/BREAKDOWNS.....
.....
.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 13 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
				Status
1. Area cleaning as per SOP.				
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.				
3. Gowning procedure is followed.				
4. Ensure that status board is displayed with mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.				
5. Ensure the absence of batch documents, labels, materials or remnants of previous product or batch on each critical part.				
6. Are the following parts of the capsule filling machine cleaned: Hopper/Feeder/Turret/Below the turret, Y-chute/Powder collection.				
7. Check and ensure that Balance Calibration & Verification records are updated.				
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40-65)% _____				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date)
				Status
1. All material from current product have been removed.				
2. All facilities have been cleaned and labeled.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Machine(s) has been checked for cleanliness.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 14 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

1. EQUIPMENT USED FOR CAPSULE FILLING

Equipment Name	Equipment ID	Mark	Equipment Name	Mark
ACG AF 90T Capsule Filling Machine	PN-37		ADU PN-38	
Capsule Mini Sorter	PN-44		EMPTY CAPS Sorter PN-41	
Capsule Polishing Machine	PN-40			
Capsule Loading Machine	PN-42			

2. MACHINE SET UP

Empty Hard Gelatin Capsule Shells (TRADE)					
Size	Colour		Printing		
# 1	Body:	Peach	Body:	KAM AMOXY 250	
	Cap:	Maroon	Cap:	KPI	
Empty Hard Gelatin Capsule Shells (UG)					
Size	Colour		Printing		
# 1	Body:	Peach	Body:	KAM AMOXY 250	
	Cap:	Maroon	Cap:	UG & KPI LOGO	
Cleaned by:		Date:	Machine speed:		
Machine setting done by:			Date:		
Checked by (section supervisor)			Date:		
Previous Product:			Batch No.:		

3. CAPSULE WEIGHT RANGE LIMITS

Average Weight of Empty Capsule Shells (mg)		% of Nominal weight	Filled Capsule wt. (mg)	Symbol
Filling Weight Per Capsule (mg)	320.00mg	Upper range	>105%	>+T2 Action
Filled Capsule Weight Range Limits (mg)			105%	=+T2 Alert
Individual Filled Capsule wt.(mg) (Nominal):			103%	+T1 Good
Accepted weight range (mg):		Norm	100%	Very Good
		Lower range	97%	-T1 Good
			95%	-T2 Alert
			<95%	<-T2 Action

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 15 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

4. CAPSULE FILLING

Started	Date:	Time:	Completed	Date:	Time:
Weight of individual filled capsules					
Average weight of empty capsule shells					

5. DETERMINE PERCENTAGE YIELD

Date:	Shift:	Operator/ Grp. Leader	Weight of filled capsules		
			Gross Wt.(kg)	Tare Wt.(kg)	Net Wt.(kg)
Total weight (kg)					

YIELD RECONCILIATION

Steps	Descriptions	Qty. (Kg)	Done By SPV (Sign & Date)	Verified By QA (Sign & Date)
A.	Weight of blend received			
B.	Actual Quantities of Capsules Filled			
C	Number Of Capsules (units)			
D	Total weight of empty shells used			
E	Sample Quantities (In-Process checks)			
F	Validation Samples			
G	Rejects (if any)			
H	Total accounted for (B+E+F)			
I	Theoretical Batch Size with respect to Blended granules = (A+D)			
Percentage Yield = $H \times 100 / I = \dots\dots\dots\%$				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):

CAUSE OF VARIATION: YES/NO

REMARKS/BREAKDOWNS.....

.....

.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 16 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
CAPSULE FILLING				

STANDARD PARAMETERS		
PARAMETERS	SPECIFICATIONS	FREQUENCIES
Appearance	Maroon cap printed 'KPI logo' in black.	Initial
	Peach body printed 'KAM AMOXY 250' in black.	Every 60 Minutes
Standard Fill		
Weight Per Capsule	$320 \pm 2\%$	Initial & Every 60 Minutes
Individual weight	$\pm 5\%$ of Average Weight	Initial
		Every 60 Minutes
Disintegration Time	Not more than 30 Minutes	Initial & Every 60 Minutes
Uniformity	$\pm 5\%$ of Average Weight. NMT 2 capsules of the 20 should be out of range of 5% of the average and none of the capsules should be out of range of 10% of the average. (pooled)	Every 60 Minutes
M/c RPM	15 rpm to 40 rpm	-----

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 17 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

6. CAPSULE FILLING IN-PROCESS QUALITY CONTROL REPORT							
M/C No.:				Location			
						Speed	
(I) Filled Capsule Weight In-Process Control							
Time HH:MM							
Done by:		Operator:		QA:			
Physical Appearance							
		Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)	Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)
Individual Weights	1						
	2						
	3						
	4						
	5						
	6						
	7						
	8						
	9						
	10						
	11						
	12						
	13						
	14						
	15						
	16						
	17						
	18						
	19						
	20						
Total wt. (mg)							
Avg wt. (mg)							
Range							
Max. wt. (mg)							
Min. Wt. (mg)							
Action >+T2							
Alert+T2							
Good +T1							
V. Good Norm							
Good -T1							
Alert -T2							
Action <-T2							
Disintegration	NMT 30min						
Lock Length(mm)	18.90-19.90						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 18 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

CAPSULE FILLING IN-PROCESS QUALITY CONTROL REPORT							
M/C No.:				Location			
						Speed	
(I) Filled Capsule Weight In-Process Control							
Time HH:MM							
Done by:		Operator:				QA:	
Physical Appearance							
		Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)	Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)
Individual Weights	1						
	2						
	3						
	4						
	5						
	6						
	7						
	8						
	9						
	10						
	11						
	12						
	13						
	14						
	15						
	16						
	17						
	18						
	19						
	20						
Total wt. (mg)							
Avg wt. (mg)							
Range							
Max. wt. (mg)							
Min. Wt. (mg)							
Action >+T2							
Alert+T2							
Good +T1							
V. Good Norm							
Good -T1							
Alert -T2							
Action <-T2							
Disintegration	NMT 30min						
Lock Length(mm)	18.90-19.90						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 19 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

CAPSULE FILLING IN-PROCESS QUALITY CONTROL REPORT							
M/C No.:				Location			
						Speed	
(I) Filled Capsule Weight In-Process Control							
Time HH:MM							
Done by:		Operator:		QA:			
Physical Appearance							
		Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)	Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)
Individual Weights	1						
	2						
	3						
	4						
	5						
	6						
	7						
	8						
	9						
	10						
	11						
	12						
	13						
	14						
	15						
	16						
	17						
	18						
	19						
	20						
Total wt. (mg)							
Avg wt. (mg)							
Range							
Max. wt. (mg)							
Min. Wt. (mg)							
Action >+T2							
Alert+T2							
Good +T1							
V. Good Norm							
Good -T1							
Alert -T2							
Action <-T2							
Disintegration	NMT 30min						
Lock Length(mm)	18.90-19.90						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 20 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

CAPSULE FILLING IN-PROCESS QUALITY CONTROL REPORT							
M/C No.:				Location			
						Speed	
(I) Filled Capsule Weight In-Process Control							
Time HH:MM							
Done by:		Operator:		QA:			
Physical Appearance							
		Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)	Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)
Individual Weights	1						
	2						
	3						
	4						
	5						
	6						
	7						
	8						
	9						
	10						
	11						
	12						
	13						
	14						
	15						
	16						
	17						
	18						
	19						
	20						
Total wt. (mg)							
Avg wt. (mg)							
Range							
Max. wt. (mg)							
Min. Wt. (mg)							
Action >+T2							
Alert+T2							
Good +T1							
V. Good Norm							
Good -T1							
Alert -T2							
Action <-T2							
Disintegration	NMT 30min						
Lock Length(mm)	18.90-19.90						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 21 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

CAPSULE FILLING IN-PROCESS QUALITY CONTROL REPORT							
M/C No.:				Location			
						Speed	
(I) Filled Capsule Weight In-Process Control							
Time HH:MM							
Done by:		Operator:		QA:			
Physical Appearance							
		Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)	Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)
Individual Weights	1						
	2						
	3						
	4						
	5						
	6						
	7						
	8						
	9						
	10						
	11						
	12						
	13						
	14						
	15						
	16						
	17						
	18						
	19						
	20						
Total wt. (mg)							
Avg wt. (mg)							
Range							
Max. wt. (mg)							
Min. Wt. (mg)							
Action >+T2							
Alert+T2							
Good +T1							
V. Good Norm							
Good -T1							
Alert -T2							
Action <-T2							
Disintegration		NMT 30min					
Lock Length(mm)		18.90-19.90					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 22 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

CAPSULE FILLING IN-PROCESS QUALITY CONTROL REPORT							
M/C No.:				Location			
						Speed	
(I) Filled Capsule Weight In-Process Control							
Time HH:MM							
Done by:		Operator:		QA:			
Physical Appearance							
		Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)	Filled Caps (mg)	Empty Shell (mg)	Net Content (mg)
Individual Weights	1						
	2						
	3						
	4						
	5						
	6						
	7						
	8						
	9						
	10						
	11						
	12						
	13						
	14						
	15						
	16						
	17						
	18						
	19						
	20						
Total wt. (mg)							
Avg wt. (mg)							
Range							
Max. wt. (mg)							
Min. Wt. (mg)							
Action >+T2							
Alert+T2							
Good +T1							
V. Good Norm							
Good -T1							
Alert -T2							
Action <-T2							
Disintegration	NMT 30min						
Lock Length(mm)	18.90-19.90						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 23 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE FILLING					

7. BULK TRANSFER / RECONCILIATION
From: Capsulation Section
To: Semi-Finished product quarantine store (SMFPQ)

Please receive the following capsules:

Date	Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)	Delivered by Operator	Received by SMFPQ Supervisor
	1					
	2					
	3					
	4					
	5					
	6					
	7					
	8					
	9					
	10					
Total						
Total number of capsules:				Average weight of capsules:		

Section supervisor (sign & date) _____

Remarks

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 24 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE INSPECTION & SORTING					

LINE CLEARANCE FORM				
Previous product:		Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date) Status
1. Area cleaning as per SOP.				
2. Ensure that all the previous product containers, material, and labels are removed from the				
3. Gowning procedure is followed.				
4. Ensure that status board is displayed with mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.				
5. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.				
6. Ensure that the Vibrosifter sieves, Metal Detector, Deduster, Scoops etc. are cleaned.				
7. Check and ensure that Balance Calibration & Verification records are updated.				
8. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____				
Remarks:				
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by QA (Sign & Date) Status
1. All material from current product have been removed.				
2. All facilities have been cleaned and labeled.				
3. The area has been cleaned and labeled.				
4. All paperwork for the current batch has been completed.				
5. Reconciliation of current batch has been done.				
6. Machine(s) has been checked for cleanliness.				
Remarks:				

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 25 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
CAPSULE INSPECTION & SORTING				

INSPECTION AND SORTING PROCEDURE

CAUTION : 1. USE SNOOD AND GLOVES WHILE HANDLING MATERIAL 2. WASH HANDS AND FACE AFTER WORKING WITH MATERIAL

1. Ensure that the inspection area and equipment are clean.
2. Ensure that all personnel are properly dressed in protective garments, gloves and masks.
3. Carry out a line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
4. Bring the capsules inside the sorting room.
5. Visually inspect the capsules for defects such as chipping, spots etc.
6. Collect all rejected capsules in properly labeled designated containers.
7. Seal inspected capsules in polythene lined drums and properly label the drums.
8. Complete the "Ending activities" section of the line clearance form.
9. Transfer inspected and sorted capsules to Packing area/Semi-finished product quarantine (SMFPQ) store.

VISUAL INSPECTION LOG

ACTIVITY	FROM	TO

RECONCILIATION AFTER INSPECTION AND SORTING:

1	Weight of capsules received for inspection and sorting _____ Kg
2	Weight of capsules received after inspection and sorting _____ Kg
3	Average weight of capsules _____ mg/capsule
4	Recyclables after inspection and sorting _____ Kg (recyclables are recorded in the BMR reconciliation sheet)
5	Total rejects of capsules after inspection and sorting _____ Kg (rejects are recorded in the BMR reconciliation sheet).
6	Total loss after inspection and sorting _____ Kg _____ %

Done by **group leader (sign & date):**

Checked by **SMFPQ supervisor (sign & date):**

Verified by **QA (sign & date):**

7	Names of person(s) involved in carrying out Inspection and Sorting		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

NOTE: Transfer inspected and sorted capsules to Semi Finished Product Quarantine (SMFPQ) Store.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 26 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
CAPSULE INSPECTION & SORTING					

QA IN-PROCESS QUALITY CONTROL REPORT FOR INSPECTION & SORTING
--

Product filled from Capsule Filling Machine Number _____

Record the net content of the drums to be sorted, weigh and record the quantity of defects sorted from the respective drums in the table below.

Drum number	Net weight (kg)	Quantity of defects (kg)	Nature of defects
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
Total			
Percentage defects = (Total quantity of defects / total net weight) x 100% = _____ %			

Done by (QA officer, sign & date): _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 27 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.		KPI/MFR/033/00		Date:	
				Shift :	

PACKAGING MATERIALS REQUISITION SHEET											
Requested by Production Supervisor			Name:		Signature:			Date:			
Approved by Production Manager			Name:		Signature:			Date:			
Pack Size: (/ mark as appropriate)					10x10 Blisters		100 Blisters X 10 Tablets		500's		1000's
S/N	Item Code	Item Description	AR. NO	Total Batch Quantity	Units	Issued by (PM Stores In-charge)		Received by (Packaging Supervisor)		Verified by (QA)	
						Sign	Date	Sign	Date	Sign	Date
FOR BLISTER PACK											
1	2AMX1CB	Kam Amoxy Capsules Aluminium Foil			Kg						
2	1AMX1CB	Kam Amoxy Capsules 10x10 Blister Cartons			Kg						
3	PF206MM	PVC Clear Film 206x0.25mm			Kg						
4	STB1XC	100X100 Capsule Blister Shipping Cartons			Kg						
5	BPTPE2	2" BOPP Printed Packing Tape			PCS						
6	PIL AMX	Kam Amoxy patient information leaflet			PCS						
7	2AMX1CB	Kam Amoxy Capsules Aluminium Foil			PCS						
8	1AMX1CB	Kam Amoxy Capsules 10x10 Blister Cartons			PCS						
9	PF206MM	PVC Clear Film 206x0.25mm			PCS						
10	STB1XC	100X100 Capsule Blister Shipping Cartons			PCS						
FOR BULK PACK											
1	3AMX4UG	Kam Amoxy capsules 1000's Labels (UG)			PCS						
2	3AMX4CJ	Kam Amoxy capsules 1000's Labels (Trade)			PCS						
3	PJ1KCC	1000 cc plastic jars with caps			PCS						
4	9LDB810	8" x 10" HDPE plastic, clear bags			PCS						
5	SJJ4XM	20 x 1000 capsules bulk shipping cartons			PCS						
6	BPTPE2	2"BOPP Printed Packing Tape			PCS						
7	PIL AMX	Kam Amoxy patient information leaflet			PCS						

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 28 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
PACKAGING					

LINE CLEARANCE FORM					
Previous product:			Batch No:		
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
1. Area cleaning as per SOP.					
2. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.					
3. Gowning procedure is followed.					
4. Ensure that status board is displayed with mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
5. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.					
6. Ensure that the Vibrosifter sieves, Metal Detector, Deduster, Scoops etc. are cleaned.					
7. Remove all the cartons, labels, strips, any material left over on the line and on the lower space of the conveyor belts.					
8. Put an authorized overprinted specimen of strip/blister/cartons/Label/Tube in the BMR for ready reference.					
9. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QAO (Sign & Date)	Status
1. All material from current product have been removed.					
2. All facilities have been cleaned and labeled.					
3. The area has been cleaned and labeled.					
4. All paperwork for the current batch has been completed.					
5. Reconciliation of current batch has been done.					
6. Stereos of previous product are submitted to production Supervisor.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 29 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
PACKAGING					

1. PACKING PROCESS & INSTRUCTIONS

BULK PACKING

1. Ensure that the packing room, packing equipment and packing table(s) are clean.
2. Carry out a Packing line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
3. Bring in the appropriate packaging material.
4. Check the labels for correctness and legibility of coding.
5. Bring in the capsules to be packed.
6. Count/Weigh the capsules.
7. Put the capsules in poly bags and seal.
8. Place the sealed poly bag in the container and firmly place the cap.
9. Re-check the batch no., mfg. date and exp. date on the labels.
10. Carefully label the container.
11. Pack the right number of containers in the appropriate corrugated shipping cartons.
12. Seal the shipping cartons with packing tape.
13. Label the shipping cartons appropriately.
14. Complete the packing records.
15. Complete the "Ending activities" section of the packing line clearance form.
16. Transfer packed goods to Finished Goods Quarantine Store and wait for QC report.

BLISTER PACKING

1. Ensure that the packing room, packing equipment and packing table(s) are clean.
2. Carry out a Packing line clearance exercise as per SOP for line clearance and complete the "beginning activities" section of the record.
3. Bring in the appropriate packaging material.
4. Set and operate the blister packing machine as per SOP for operating blister packing machine.
5. Check the first code for correctness and legibility of coding.
6. Bring in the tablets to be blistered and complete the operation and transfer all those blisters to secondary packing area.
7. Pack 10 blisters per unit carton/ 100 blisters per unit dispenser.
8. Pack the right number of unit cartons in the corrugated shipping cartons.
9. Seal the shipping carton with packing tape.
10. Label the shipping carton appropriately.
11. Complete the packing records.
12. Complete the "Ending activities" section of the packing line clearance form.
13. Transfer packed goods to Finished Goods Quarantine Store and wait for QC report.

Reference Documents:

1. SOP for status labelling.
2. SOP for line clearance.
3. SOP for Bulk packing.
4. SOP for Blister packing

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 30 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
PACKAGING					

2. CAPSULE BLISTERING

1. BLISTER MACHINE SETTING

Blister Machine number		Blister machine Make	
PVC Clear Film Size		Aluminium Foil Size	
Set Sealing Temperature		Observed Sealing Temperature	
Set Forming Temperature		Observed Forming Temperature	
Machine Speed		Blister Pack Size	
Machine Setting Done By (Sign & Date):			
Checked by Supervisor (Sign & Date):			
Verified by QA (Sign & Date):			

2. BLISTER PACKING

Started	Date:	Time:	Completed	Date:	Time:
---------	-------	-------	-----------	-------	-------

3. DETERMINE PERCENTAGE YIELD

Date:	Shift:	Operator/ Grp. Leader	Number of Blisters
Total Number of Blisters			

YIELD RECONCILIATION

Steps	Descriptions	Qty. (units)	Checked By SPV (Sign & Date)	Verified By QA (Sign & Date)
A.	Theoretical Batch Size (Blisters) with respect to capsules received from sorting.			
B.	Actual Quantities of Blisters			
C.	Sample Quantities (In-Process checks)			
D.	Validation Samples			
E.	Rejects (if any)			
F.	Total Actual Yield (B+C+D)			
G.	Unaccountable Losses [A-(E+F)]			
Percentage Yield = $F \times 100 / A = \dots\dots\dots\%$				
Permissible Yield Limit: (98% to 102%)				

(If percentage yield is within 98% to 102% no justification is required, otherwise explain variation in space below):

CAUSE OF VARIATION: YES/NO

REMARKS/BREAKDOWNS.....

.....

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 31 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
PACKAGING					

3. PACKING CODING CONTROL

1	Batch No:	Coding approved by Packing Supervisor (sign & date)	
2	Mfg. Date:		
3	Exp.Date:		
4	Item(s) to be coded: Labels, Unit Cartons, Shipper labels, Blisters, Others (specify); mark as appropriate		
5	Coding Stamp set by (operator): sign & date:		
	Coding Machine Number:		
6	Attach first coded samples and patient information leaflet (PIL) here:		
7	Coding approved by: QA Officer (Sign & Date):		
8	Names of person(s) applying labels and packing - BULK & BLISTER (mark as appropriate).		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)
9	Coding Reconciliation:		
	Name of item:		
	No. of item issued:		
	No. of item coded:		
	No. of item uncoded:		
	No. of items damaged / rejected during coding:		
	No. of items used:		
	No. of excess items, if any:		
	No. of item to be destroyed (with QA consent):		
	Reconciliation approved by: QA (Sign & Date):		

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 32 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
PACKAGING					

4. BULK TRANSFER AFTER INSPECTION & SORTING
--

From: Semi finished product quarantine store (SMFPQ)

To: Packing Department (Bulk / blisters/strips)

Please receive the following capsules for packing:

Date	Drum No.	Gross Weight (Kg)	Tare Weight (Kg)	Net Weight (Kg)	Delivered by SMFPQ supervisor	Received by Operator/ Packing supervisor
	1					
	2					
	3					
	4					
	5					
	6					
	7					
	8					
	9					
	10					
Total						
Total number of capsules:				Average weight of capsules:		

Quantity packed (bulk /blisters): -----

(A) REJECTS AFTER PACKING _____ KG/G

(B) RE-WORKABLES AFTER PACKING _____ KG/G

Section supervisor (sign & date) _____

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 36 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
SECONDARY PACKAGING					

LINE CLEARANCE FORM					
Previous product:		Batch No:			
BEGINNING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements.		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	
				QA (Sign & Date)	Status
10. Area cleaning as per SOP.					
11. Ensure that all the previous product containers, material, and labels are removed from the manufacturing area.					
12. Gowning procedure is followed.					
13. Ensure that status board is displayed with mentioning Product Name, Batch No, B. Size, Mfg. Date, Exp.Date & Status with sign & date.					
14. Check the proper cleanliness of the Air Supply and Return grills and Air Conditioning System.					
15. Ensure that the Vibrosifter sieves, Metal Detector, Deduster, Scoops etc. are cleaned.					
16. Remove all the cartons, labels, strips, any material left over on the line and on the lower space of the conveyor belts.					
17. Put an authorized overprinted specimen of strip/blister/cartons/Label/Tube in the BMR for ready reference.					
18. Ensure Environmental conditions are met as per SOP. Temperature NMT 28 °C _____ Relative Humidity (40 to 65%) _____					
Remarks:					
ENDING ACTIVITIES: The following must be checked to ensure adherence to GMP requirements		Done by Operator (Sign & Date)	Checked by Supervisor (Sign & Date)	Verified by	Status
				QAO (Sign & Date)	
7. All material from current product have been removed.					
8. All facilities have been cleaned and labeled.					
9. The area has been cleaned and labeled.					
10. All paperwork for the current batch has been completed.					
11. Reconciliation of current batch has been done.					
12. Stereos of previous product are submitted to production Supervisor.					
Remarks:					

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 37 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
SECONDARY PACKAGING				

INSPECTION AND SORTING OF BLISTERS AND STRIPS
--

INSPECTION AND SORTING PROCEDURE

1. Ensure Inspection and sorting area is clean as per the procedure for Line Clearance ,SOP No. GEN/013.
2. Carry out status labelling of the area and facilities if any as per the procedure for Status Labelling and control of Status Labels, SOP NO. GEN/002.
3. Bring the strips/blisters to the prepared Inspection and Sorting area in the Secondary Packing section.
4. Visually inspect/Check the strips / blisters / coded materials for all the possible defects. Some examples of types of rejection on packing line are Incorrect coding, Text matter/hologram/smudging of ink over printed foil, Rejected strips due to cutting, Empty pocket, Pocket puncture, Sealing rejection, Overheated strips/blisters, Chipped / capped / broken tablets, Improper forming of blisters, Improper appearance / absence of hologram in the final blister strip among others.
5. After the Inspection and Sorting operation, the Supervisor shall intimate to the QA Officer online to verify the integrity of the sorted blisters/strips and approve packing to start.
6. Remove all defective blisters/strips, place them in a polybag or drum and label as recyclables. These should be de-blistered/de-stripped once sorting and inspection of blisters is done.

NOTE:

1. All personnel involved in Inspection and Sorting shall undergo visual checks (Ocular test) annually to ascertain proper vision with or without glasses and Color blindness.
2. Report to the Packing / Sorting Supervisor or Production Manager in case you start to feel uneasy or have a headache during the inspection and sorting process.
3. Personnel conducting the inspection and sorting shall where applicable take a break every after 3 hours to refresh their eyesight.

Inspection and Sorting Line Clearance Done by (sign & date)	
Checked by Production Supervisor (sign & date):	
Verified by QA (sign & date):	

VISUAL INSPECTION LOG		
ACTIVITY	FROM	TO

10	Names of person(s) involved in carrying out Inspection and Sorting		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 38 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
SECONDARY PACKAGING					

SECONDARY PACKING PROCEDURE			
1. Check the specimen attached to the BMR of overprinted carton/Monocartons, leaflet, and shipper label if it has been checked by the line supervisor and approved by the QA Officer online. 2. Transfer previously approved overprinted packing material from the printing room to the secondary packing room. 3. Check if displayed details on the carton/Monocartons/dispensers of the batch number, manufacturing and expiry date are as per BMR/BPR. Check carton for overprinting intactness. 4. Check if the information on leaflets/inserts is clear, legible and specific to the product to be packed. 5. Open the carton and insert leaflet with number of strips/blisters as specified in the BMR. Close the carton and arrange it into the shipper. 6. Once the shipper is filled with the number of cartons as specified in the BMR, insert a fully filled packing slip with name of the product & batch number, the name of the packing personnel and this should be checked and signed by the group leader. 7. Group leaders should check the shipper to ensure that it has the correct quantity of cartons/Monocartons/dispensers. If this is okay, group leader shall sign the packing slip and intimate to the Supervisor and QA Officer online to confirm the correctness of the quantity in the shipper and other packing requirements. The Supervisor and the QA officer shall then sign on the packing slip and authorize the group leader to seal the shipper with BOPP tape and attach the shipper label. Exact quantity of loose shipper is written on shipper label and sealed. 8. Once the entire batch has been completed, transfer to finished product quarantine store. 9. Any leftover printed cartons/Monocartons/dispensers/labels should be destroyed in presence of QA Officer online and the number recorded in the BMR. 10. Transfer leftover packaging material (that is not printed on) to Packaging Material Stores with a return note.			
Secondary Packing Line Clearance Done by (sign & date)			
Checked by Production Supervisor (sign & date):			
Verified by QA (sign & date):			
RECONCILIATION AFTER PACKING			
Product Name:		Batch Number:	
Batch Size (Cartons/Jars):		Pack Size:	
No. of intact shippers packed		No. of cartons/jars in intact Shippers	
Loose Quantity Packed:		Total Quantity Packed (cartons/jars)	
Done by group leader (sign & date):			
Checked by SMFPQ supervisor (sign & date):			
Verified by QA (sign & date):			
11	Names of person(s) applying labels and packing - BULK & BLISTER (mark as appropriate).		
	(a)	(b)	(c)
	(d)	(e)	(f)
	(g)	(h)	(i)

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY		Page 39 of 45
BATCH SIZE		500,000 Capsules	Mfg. Date :	
BATCH No.			Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:	Shift :	
SECONDARY PACKAGING				

SHIPPER WEIGHT VERIFICATION

- After the packing operation is completed for the batch, **Weigh the packed shippers one by one** and record the weight of packed shipper in BPR and also on the Shipper Label and recheck the shipper if any discrepancy is found.

Average Weight of Empty Carton:					Average Weight of Packed Unit Carton:						
Average Weight of Empty Shipper:					Average Weight of Packed Shipper:						
Acceptance Limit: $\pm 2.0\%$ of average weight of packed Shipper					Range: _____ to _____						
SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT	SHIPPER No.	SHIPPER WEIGHT
01		11		21		31		41		51	
02		12		22		32		42		52	
03		13		23		33		43		53	
04		14		24		34		44		54	
05		15		25		35		45		55	
06		16		26		36		46		56	
07		17		27		37		47		57	
08		18		28		38		48		58	
09		19		29		39		49		59	
10		20		30		40		50		60	

Remarks: All shippers are within Specifications YES/NO _____

Weighing Done by (Sign & Date):

Checked by Packing Supervisor (Sign & Date):

Approved by QA (Sign & Date):

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 43 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Shift :	
FINISHED PRODUCT RECONCILIATION SHEET					

Quantity Packed	Retention Sample Quantity	Actual Quantity Released	Quantity Delivered to FGS	Delivery Note Number	Delivery Date

FINAL BATCH RECONCILIATION

ITEM	Quantity (Gm/Kg).	Quantity* (No. of Units)	Percent =[QTYx100/(A)]
(A) Theoretical Yield =(Batch Size):			
(B) Yield - 1st Crop (Quantity Packed):			
(C) In-Process/QC Samples:			
(D) Re-workables =[(I)+(II)]:			
(I) Powder/Granules*			
(II) Capsules			
(E) Rejects =[(I)+(II)]:			
(I) Powder/Granules*			
(II) Capsules			
(F) Total Yield =(B+C+D):			
(H) Shortage =[A-(E+F)]:			

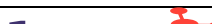
Done by Supervisor (Sign & Date):	Verified by QAO (Stamp, Sign & Date):
--	--

Reconciliation Approved By **Production Manager (Sign & Date):** _____

*Quantity (No. of units) is derived by converting this quantity using average fill weight.

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.	
PRODUCT		KAM AMOXY		Page 44 of 45	
BATCH SIZE		500,000 Capsules		Mfg. Date :	
BATCH No.				Exp. Date :	
MFR No.	KPI/MFR/033/00	Date:		Actual Quantity :	
BATCH RELEASE SHEET					

The following Batch production records have been properly filled out, checked, dated and signed by authorized person(s). The tests are also within limits. The batch is hereby approved for release.		
Reviewed and approved by QA Head/Designee (Sign & Date):		
Records Reviewed (Mark as appropriate)	QA Head/Designee	Company /Prod. Pharmacist
1. Raw Material Requisition/Dispensing Sheet		
2. Mixing Records		
3. Blending / Lubrication Records		
4. Capsule filling/ In-process Control Records		
5. Process Line Clearance Records		
6. Records for Semi-finished Products (blended material, filled capsules etc.)		
7. Inspection and Sorting Records		
8. Packing Line Clearance Records		
9. Packing Material Authorization Note		
10. Packing Coding Control Records		
11. Product Reconciliation Sheet		
12. IPQA Reports		
13. QC Analytical reports		
13.1. Finished Product Analytical Report		
13.2. Finished Product Analytical Data Sheet (Observations, Calculations, Printouts, References).		
I hereby certify that all the manufacturing stages of this batch of finished product have been carried out in full compliance with the Uganda National Drug Authority (NDA) GMP requirements and with the requirements of the product registration. I, therefore, AUTHORISE release of the batch.		
Company Pharmacist/Pharmacist In-Charge:		
_____	_____	_____
Name	Signature	Date

		Batch Manufacturing Record and Procedure		KAMPALA PHARMACEUTICAL INDUSTRIES (1996) LTD.
PRODUCT		KAM AMOXY	Page 45 of 45	
BATCH SIZE		500,000 Capsules		
MFR No.	KPI/MFR/033/00			
REVISION HISTORY				

REVISION NO.	CHANGES INCORPORATED	REASON OF CHANGE	EFFECTIVE DATE
01	- KPI new logo to be incorporated - The provision for BMR preparation, Review, Approval and Authorization have been introduced on the Master Formula Page(s). - An Index has been introduced. - The Abbreviations used in the BMR have been defined in full. - List Of Equipment & Instruments used in product manufactured have been summarized. - Line Clearance for Dispensing previously missing has been introduced. - Documentation of First Drying, Second Drying, and Final Drying have been split and captured independently, and Milling activities including Sieves used plus LOD data are also being documented at granulation. Identity of the FBD used to dry a particular Lot has been captured. - In-process Quality Control and Assurance (IPQC/IPQA) reports for both Production and QA at all stages (as appropriate) have been incorporated. - Initial Compression Machine Setting instructions, Initial Dies and Punches Checking Record, Initial Individual Weight Records as per number of Punches/Stations of the Compression machine have been introduced/Captured.	To improve on clarity of the BMR. To match up with the current practices.	
02	- To eliminate duplication in the BMR header, introduce page breaks for different process steps and include MFR numbers in the BMR. - To include the Line Clearance for Inspection and Sorting of Blisters/Strips and the Secondary Packing Procedure.	To improve coherence of the BMR To ensure maximum elimination of blister/Strip defects before packing.	
03	- Batch Release sheet revised to ensure that batch release is Authorized by the Company Pharmacist/Pharmacist In-Charge. - Shipper weight verification record introduced. - Packaging Material Requisition Sheet incorporated in the BMR.	To comply with regulatory requirements for batch release. To eliminate the risk of having less or more product quantity in a shipper	